

Innovative Biomaterials for Retinal Surgery

Ultrapure fluids and gases
made in Germany



Content

Solutions for ophthalmic surgery



G-81010
 G-81005
 G-81015
 G-81025
 G-81002
 G-81001



5 ml: G-80315
 7 ml: G-80317
 5 ml: G-80305
 7 ml: G-80307
 5 ml: G-80115
 7 ml: G-80117
 5 ml: G-80105
 7 ml: G-80107



G-80950
 G-80960
 G-80970

Dyes

Posterior segment

Brilliant Peel®	4
Brilliant Peel® Dual Dye	8

Anterior segment

Vioron®	12
---------------	----



Intraoperative Tamponades


F-Decalin.....	14
F-Octane	14




Long-Term Tamponades

Gases


EasyGas® SF ₆	16
EasyGas® C ₃ F ₈	16
EasyGas® C ₂ F ₆	16

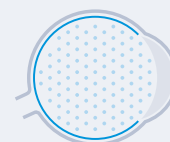

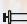

 G-80740
  G-80750
  G-80720
  G-80820
 G-80710
  G-80810

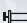

 G-80925



 G-80615
  6 ml: G-80606
 9 ml: G-80609

Long-Term Tamponades

Silicone oils


Siluron® Xtra	18
Siluron® 2000	18
Siluron® 1000	19
Siluron® 5000	19



Long-Term Tamponades

Heavy silicone oils

Densiron® Xtra	22
----------------------	----

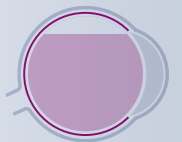

Cleansing Fluids / WashOut

F4H5® WashOut.....	26
F6H8® Vitreous Substitute	27



Accessories

Accessoires PFCL / EasyGas®	28
Accessories silicone oils.....	29

Fluoron

Purity and variety made in Germany... 31

Brilliant Peel®

Heavy dye for selective staining
of the ILM



Precise and intense staining
of the ILM due to fast sinking dye

Safe application under air and BSS

Quick and easy application
(Ready-to-use)

Physiological osmolarity

Biocompatible

Composition and properties

1 ml of Brilliant Peel® contains:

- 0.25 mg Brilliant Blue G
- Disodium phosphate ($\text{Na}_2\text{HPO}_4 \times 2 \text{H}_2\text{O}$)
- Monosodium phosphate ($\text{NaH}_2\text{PO}_4 \times 2 \text{H}_2\text{O}$)
- Sodium chloride (NaCl)
- Deuterium oxide (D_2O)
- Water for injection

Density: 1.02 g/cm³

Packaging units




G-81010 Brilliant Peel® Syringe
0.5 ml Syringe, 5 pcs. per box, sterile



G-81005 Brilliant Peel® Vial
0.5 ml Vial, 5 pcs. per box, sterile

Fields of application

Brilliant Peel® was developed for specific staining of the inner limiting membrane (ILM). Specific staining of the ILM allows it to be clearly differentiated from the underlying retinal tissue and the epiretinal membranes. Due to the density of 1.02 g/cm³ Brilliant Peel® quickly sinks to the fundus of the eye without diffuse dispersion in the whole bulbus. The surgically demanding removal of the ILM thus becomes easier and safer.

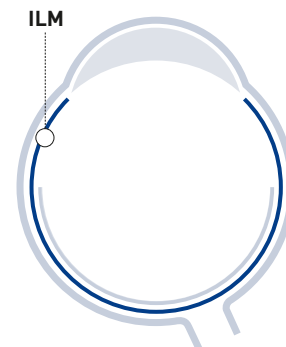

Comparison of Brilliant Blue G (BBG), Indocyanine Green (ICG) and Trypan Blue (TB) for chromovitrectomy^{22, 25}

	BBG	ICG	TB
Chemical classification	Triphenylmethane	Cyanine	Diazo
Color	blue	dark green	dark blue
Ready-to-use	yes	no	yes
Toxicity ^{12, 16, 17, 20, 21, 23}	no	yes	moderate
Significant cytotoxic effect ^{23, 30} (according to DIN EN ISO 10993)	> 0.3 g/l reduces cell growth	> 0.24 g/l causes apoptosis	> 0.13 g/l
Approval	yes	no	yes
Affinity for ILM ^{18, 19, 22}	high	high	low
Affinity for ERM ¹²	low	low	high
Selective staining of ILM ^{12, 18, 19}	high	high	low
Exposure time	short	short	long
Liquid / Gas exchange necessary	no	partially	yes

Testimonials Brilliant Peel®

”

“Our data underline the good biocompatibility of BBG and its applicability and safety for the use in humans. BBG provides a sufficient and selective staining of the ILM. No retinal toxicity related to BBG was observed in our animal study and our shortterm clinical investigation in humans.”

Remy, M., S. Thaler, R. G. Schumann et al. 2008. "An in vivo evaluation of Brilliant Blue G in animals and humans" British Journal of Ophthalmology 92(8): 1142-1147.

”

“Heavy brilliant blue G (BBG-D₂O) provides a significantly improved staining effect of the ILM and by this makes ILM peeling more efficient, easier, faster and less traumatic.”

Gerding, H., M. Timmermann and U. Thelen. 2011. "Intravital staining of the internal limiting membrane with a novel heavy solution of brilliant blue G." Klinische Monatsblätter für Augenheilkunde, 228(04): 298-301.

”

“Brilliant blue G-D₂O dye compartment is convenient, as the dye sinks readily onto the retinal surface and dye dispersion to the remaining vitreous is reduced. Indications for dye-related toxicity or complications were not seen.”

Henrich, P. B., C. Valmaggia, C. Lang, S. G. Priglinger, C. Haritoglou, R. W. Strauss and P. C. Cattin. 2013. "Contrast recognizability during Brilliant Blue G - and heavier-than-water Brilliant Blue G-assisted chromovitrectomy: a quantitative analysis." Acta Ophthalmologica 91(2): e120-124.

”

“Although the MH closure rate was the same using BBG or ICG for ILM peeling, visual acuity improvement was better in eyes peeled with BBG compared to eyes peeled with ICG.”

Jenisch, T. M., F. Zeman, M. Koller, D. A. Märker, H. Helbig and W. A. Herrmann. 2017. "Macular hole surgery: An analysis of risk factors for the anatomical and functional outcomes with a special emphasis on the experience of the surgeon." Clinical Ophthalmology (Auckland NZ) 11: 1127-1134.

References

1. Remy M., S. Thaler, R. G. Schumann, C. A. May, M. Fiedorowicz, F. Schuettauf, M. Grünerich, S. G. Priglinger, M. M. Nentwich, A. Kampik et al. 2008. "An in vivo evaluation of Brilliant Blue G in animals and humans" *British Journal of Ophthalmology* 92(8): 1142-1147.
2. Henrich, P. B., C. Valmaggia, C. Lang and P. C. Cattin. 2014. "The price for reduced light toxicity: Do endoillumination spectral filters decrease color contrast during Brilliant Blue G-assisted chromovitrectomy?" *Graefe's Archive for Clinical and Experimental Ophthalmology* 252(3): 367-74.
3. Henrich, P. B., C. Valmaggia, C. Lang, S. G. Priglinger, C. Haritoglou, R. W. Straus and P. C. Cattin. 2013. "Contrastrecognisability during Brilliant Blue G - and heavier-than-water Brilliant Blue G-assisted chromovitrectomy: a quantitative analysis." *Acta Ophthalmologica* 91(2): e120-4.
4. Williamson, T. H. and E. Lee. 2014. "Idiopathic macular hole: analysis of visual outcomes and the use of Indocyaninegreen or Brilliant Blue for internal limiting membrane peel." *Graefe's Monatsblätter für Augenheilkunde* 230(7): 701-6.
5. Ooi Y. L., T. F. Khang, M. Naidu and K. C. Fong. 2013. "The structural effect of intravitreal Brilliant Blue G and Indocyanine Green in rats eyes." *Eye (London)* 27(3): 425-31.
6. Baba, T., A. Hagiwara, E. Sato, M. Arai, T. Oshitari and S. Yamamoto. 2012. "Comparison of vitrectomy with BrilliantBlue G or Indocyanine Green on retinal microstructure and function of eyes with macular hole." *Ophthalmology* 119(12): 2609-2615.
7. Januschowski, K., S. Mueller, M. S. Spitzer, C. Schramm, D. Doycheva, K. U. Bartz-Schmidt and P. Szurman "Evaluating retinal toxicity of a new heavy intraocular dye, using a model of perfused and isolated retinal cultures of bovine and human origin." *Graefe's Archive for Clinical Experimental Ophthalmology* 250(7): 1013-1022.
8. Pelayes, D. E., F. Kuhn, A. M. Folgar, W. Takahashi, A. Bastien, P. N. Vinicius and J. O. Zarate. 2012. "Staining of the internal limiting membrane with the use of heavy Brilliant Blue G." *Ophthalmic Research*, 48 [Suppl. 1]: 21-25.
9. Gerding H., M. Timmermann and U. Thelen. 2011. "Intravitreal staining of the internal limiting membrane with a novel heavy solution of Brilliant Blue G." *Klinische Monatsblätter für Augenheilkunde* 228(4): 298-301.
10. Haritoglou, C., R. G. Schumann, A. Kampik and A. Gandorfer "Heavy Brilliant Blue G for internal limiting membrane staining." *Retina* 31(2): 405-7.
11. Haritoglou, C. and A. Kamik. 2006. "Färbetechniken in der Makulachirurgie." *Der Ophthalmologe* 103: 927-934.
12. Höing, A., M. Remy, M. Dirisamer, S. Priglinger, C. L. Schönfeld, A. Kampik and C. Haritoglou. 2011. "An in-vivo evaluation of Brilliant Blue G in macular surgery." *Klinische Monatsblätter für Augenheilkunde* 228(8): 724-728.
13. Jenisch, T. M., F. Zeman, M. Koller, D. A. Märker, H. Helbig and W. A. Herrmann. 2017. "Macular hole surgery: An analysis of risk factors for the anatomical and functional outcomes with a special emphasis on the experience of the surgeon." *Clinical Ophthalmology [Auckland NZ]* 11: 1127-1134.
14. Wirbelauer, C., H. Häberle and D. T. Pham. 2011. "Clinical experience with Brilliant Blue G staining of the retinal surface." *Klinische Monatsblätter für Augenheilkunde* 228(1): 62-65.
15. Lüke, C., M. Lüke, T. S. Dietlein, A. Hueber, J. Jordan, W. Sickel and B. Kirchhof. 2005. "Retinal tolerance to dyes." *British Journal of Ophthalmology* 89(9): 1188-1191.
16. Ueno, A., T. Hisatomi, H. Enaida, T. Kagimoto, Y. Mochizuki, Y. Goto, T. Kubota, Y. Hata and T. Ishibashi. 2007. "Biocompatibility of Brilliant Blue G in a rat model of subretinal injection" *Retina* 27(4): 499-504.
17. Enaida, H., T. Hisatomi, Y. Hata, A. Ueno, Y. Goto, T. Yamada, T. Kubota and T. Ishibashi. 2006. "Brilliant Blue G selectively stains the internal limiting membrane - Brilliant Blue G assisted membrane peeling." *Retina* 26(6): 631-636.
18. Enaida, H., T. Hisatomi, Y. Goto, Y. Hata, A. Ueno, M. Miura, T. Kubota and T. Ishibashi. 2006. "Preclinical investigation of internal limiting membrane staining and peeling using intra-vitreous Brilliant Blue G." *Retina* 26(6): 623-630.
19. Hisatomi, T., H. Enaida, H. Matsumoto, T. Kagimoto, A. Ueno, Y. Hata, T. Kubota, Y. Goto and T. Ishibashi. "Staining ability and biocompatibility of Brilliant Blue G - preclinical study of Brilliant Blue G as an adjunct for capsular staining." *Archives of Ophthalmology* 124(4): 514-519.
20. Goldman, J. M. and C. L. Karp. 2007. "Adjunct devices for managing challenging cases in cataract surgery - capsular staining and ophthalmic viscosurgical devices." *Current Opinion in Ophthalmology* 18(1): 52-57.
21. Meyer, C. H., R. A. Oechsler and E. B. Rodrigues. 2007. "Historical considerations in applying vital dyes in vitreoretinal surgery: from early experiments to advanced chromovitrectomy." *Expert Review of Ophthalmology* 2(1): 71-77.
22. Hiebl, W., B. Günther and H. Meinert. 2005. "Substances for staining biological tissues: use of dyes in ophthalmology." *Klinische Monatsblätter für Augenheilkunde* 222(4): 309-311.
23. Schuettauf, F., C. Haritoglou, Ch. A. May, R. Rejdak, A. Mankowska, W. Freyer, K. Eibl, E. Zrenner, A. Kampik, and S. Thaler. 2006. "Administration of Novel Dyes for Intra-ocular Surgery: An In Vivo Toxicity Animal Study." *Investigative Ophthalmology & Visual Science* 47(8): 3573-3578.
24. Rodrigues, E. B., M. Maia, C. H. Meyer, F. M. Penha, E. Dib and M. E. Farah. 2007. "Vital dyes for chromovitrectomy." *Current Opinion in Ophthalmology* 18(3): 179-87.
25. Furlani, B. A., L. Barroso, D. Sousa-Martins, M. Maia, M. N. Moraes-Filho, E. Badaro, R. Portella, A. A. Lima-Filho, E. Buchelle Rodrigues and R. Belfort. 2014. "Lutein and zeaxanthin toxicity with and without brilliant blue in rabbits." *Journal of Ocular Pharmacology and Therapeutics* 30(7): 559-566.
26. Tura, A., A. Alt, C. Haritoglou, C. H. Meyer, T. Schneider, S. Grisanti, J. Lüke and Matthias Lüke. 2014. "Testing the effects of the dye acid violet-17 on retinal function for an intraocular application in vitreo-retinal surgery." *Graefe's Archive for Clinical and Experimental Ophthalmology* 252(12): 1927-1937.
27. Sousa-Martins, D., M. Maia, M. Moraes, A. A. Lima-Filho, E. B. Rodrigues, J. Chen, M. E. Farah, L. B. Santos and R. Belfort Jr. 2012. "Use of lutein and zeaxanthin alone or combined with Brilliant Blue to identify intra-ocular structures intraoperatively." *Retina*. 32(7): 1328-1336.
28. Patent DE102012103097 A1 30. Kawahara S, et al. 2007. "Intracellular events in retinal glial cells exposed to ICG and BBG." *IOVS. Vol. 48, No. 10*

Brilliant Peel® Dual Dye

The non-toxic dual dye



Safer peeling due to distinct staining of the membrane



Intense and selective staining of ILM, ERM and vitreous remnants

Fast sinking – maximized contact surface with tissue due to higher density

Safe application under air and BSS

Quick and easy application (Ready to use)

Physiological osmolarity

Biocompatible


Composition and properties


1 ml of Brilliant Peel® Dual Dye contains:

- 0.25 mg Brilliant Blue G
- 1.3 mg Bromphenol Blue
- Disodium phosphate ($\text{Na}_2\text{HPO}_4 \times 2 \text{H}_2\text{O}$)
- Monosodium phosphate ($\text{NaH}_2\text{PO}_4 \times 2 \text{H}_2\text{O}$)
- Sodium chloride (NaCl)
- Deuterium oxide (D_2O)
- Water for injection

Density: 1.03 g/cm³

Packaging units

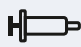 **G-81015 Brilliant Peel® Dual Dye Syringe**
0.5 ml syringe, 5 pcs. per box, sterile

 **G-81025 Brilliant Peel® Dual Dye Vial**
0.5 ml vial, 5 pcs. per box, sterile


Video

Scan QR-Code for further information on Brilliant Peel Dual Dye



Fields of application

Brilliant Peel® Dual Dye was developed for specific staining of the inner limiting membrane (ILM) and epiretinal membrane (ERM). Specific staining of the ILM and ERM allows them to be clearly distinguished from the underlying retinal tissue, thus making the challenging surgical removal of the ILM and ERM easier and safer. Due to the density of 1.03 g/cm³ Brilliant Peel® Dual Dye quickly sinks to the fundus of the eye without diffuse dispersion in the whole bulbus. Brilliant Peel® Dual Dye can, to some extent, also be used for staining vitreous remnants.

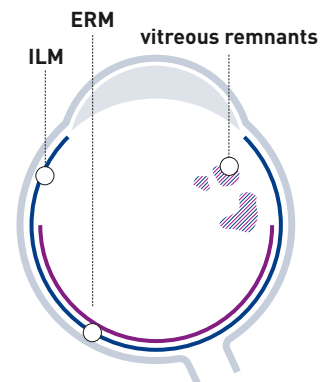

Comparison of Brilliant Blue G (BBG), Bromphenol Blue (BPB), Indocyanine Green (ICG), Trypan Blue (TB) and Lutein for chromovitrectomy^{38, 39}

	Brilliant Peel® Dual Dye	Other dyes		
	BBG & BPB	ICG	TB	Lutein
Chemical classification	Triphenylmethane	Cyanine	Diazo	Carotenoid
Color	violet-blue	green	blue	yellow-orange
Dyes⁴²	Brilliant Blue G & Bromphenol Blue	Indocyanine Green	Trypan Blue	Lutein
Toxicity^{31, 32, 33, 36, 37, 39, 40, 42, 45}	no	yes	moderate	no
Approval	yes	no	yes	yes
Affinity for ILM^{34, 35, 38, 43, 44, 45}	high	high	low	no
Affinity for ERM^{32, 43, 45}	high	low	high	n.a.
Affinity for vitreous remnants^{41, 43, 44, 45}	moderate	low	low	high
Exposure time	short	short	long	short
Liquid / Gas exchange	no	no	yes	no

Testimonials Brilliant Peel® Dual Dye



”

”Excellent staining of pre-retinal membranes and vitreous remnants.“

Senior Consultant Dr. Jürgen Steinhauer, Eye Clinic of Helios University Hospital Wuppertal – University Witten/Herdecke, Germany



”

”Outstanding staining properties and an impressive sinking behavior makes Brilliant Peel Dual Dye the perfect tool for a safe peeling in epiretinal macular procedures. Flawless for a fast and reliable multiple staining of different membrane parts.“

Prof. Dr. Lars-Olof Hattenbach, Director of Eye Clinic Ludwigshafen, Germany



”

”Even under yellow UV-IOL the shape of the retinal nerve fiber layer (RNFL) on the ILM was perfectly visible. A highly promising dye with excellent sinking properties.“

Senior Physician Dr. A. Viestenz, University Clinic of Saarland, Homburg, Germany

References

31. Lüke, C., M. Lüke, T. S. Dietlein, A. Hueber, J. Jordan, W. Sickel and B. Kirchhoff. 2005. "Retinal tolerance to dyes". *British Journal of Ophthalmology*, 89 (9): 1188-1191.
32. Haritoglou, C. and A. Kampik. 2006. "Färbetechniken in der Makulachirurgie". *Der Ophthalmologe* 103(11): 927-934.
33. Ueno, A., T. Hisatomi, H. Enaida, T. Kagimoto, Y. Mochizuki, Y. Goto, T. Kubota, Y. Hata and T. Ishibashi. 2007. "Biocompatibility of Brilliant Blue G in a rat model of subretinal injection". *Retina* 27(4): 499-504.
34. Enaida, H., T. Hisatomi, Y. Hata, A. Ueno, Y. Goto, T. Yamada, T. Kubota and T. Ishibashi. 2006. "Brilliant Blue G selectively stains the internal limiting membrane - Brilliant Blue G assisted membrane peeling." *Retina* 26(6): 631-636.
35. Enaida, H., T. Hisatomi, Y. Goto, Y. Hata, A. Ueno, M. Miura, T. Kubota, T. Ishibashi. 2006. "Preclinical investigation of internal limiting membrane staining and peeling using intravitreal Brilliant Blue G." *Retina* 26(6): 623-630.
36. Hisatomi, T. H. Enaida, H. Matsumoto, T. Kagimoto, A. Ueno, Y. Hata, T. Kubota, Y. Goto, T. Ishibashi. 2006. "Staining ability and biocompatibility of Brilliant Blue G - preclinical study of Brilliant Blue G as an adjunct for capsular staining." *Archives of Ophthalmology* 124(4): 514-519.
37. Goldman, J. M. and C. L. Karp. 2007. "Adjunct devices for managing challenging cases in cataract surgery - capsular staining and ophthalmic viscosurgical devices." *Current Opinion in Ophthalmology* 18(1): 52-57.
38. Meyer, C. H., A. Oechster and E. B. Rodrigues. 2007. "Historical considerations in applying vital dyes in vitreoretinal surgery: from early experiments to advanced chromovitrectomy." *Expert Review of Ophthalmology* 2(1): 71-77.
39. Hiebl, W., B. Günther and H. Meinert. 2005. "Substances for staining biological tissues: use of dyes in ophthalmology." *Klinische Monatsblätter für Augenheilkunde* 222(4): 309-311.
40. Schuettauf, F., Ch. Haritoglou, Ch. A. May, R. Rejdak, A. Mankowska, W. Freyer, K. Eibl, E. Zrenner, A. Kampik and S. Thaler. 2006. "Administration of Novel Dyes for Intraocular Surgery: An In Vivo Toxicity Animal Study." *Investigative Ophthalmology & Visual Science* 47(8): 3573-3578.
41. Rodrigues, E. B., M. Maia, C. H. Meyer, F. M. Penha, E. Dib and M. E. Farah. 2007. "Vital dyes for chromovitrectomy." *Current Opinion in Ophthalmology* 18(3): 179-87.
42. Furlani, B. A., L. Barroso, D. Sousa-Martins, M. Maia, M. N. Moraes-Filho, E. Badaro, R. Portella, A. A. Lima-Filho, E. Buchelle Rodrigues and Rubens Belfort Jr. 2014. "Lutein and zeaxanthin toxicity with and without brilliant blue in rabbits." *Journal of Ocular Pharmacology and Therapeutics* 30(7): 559-566.
43. Patent DE102012103097 A1
44. Sousa-Martins, D. M. Maia, M. Moraes, A. A. Lima-Filho, E. B. Rodrigues, J. Chen, M. E. Farah, L. Santos and Rubens Belfort Jr. 2012. "Use of lutein and zeaxanthin alone or combined with Brilliant Blue to identify intraocular structures intraoperatively." *Retina* 32(7): 1328-1336.
45. Haritoglou, C., R. G. Schumann, R. Strauss, S. G. Priglinger, A. S. Neubauer and A. Kampik. 2007. "Vitreoretinal surgery using bromphenol blue as a vital stain: evaluation of staining characteristics in humans." *British Journal of Ophthalmology* 91(9): 1125-1128.
46. Kodjikian, L., T. Richter, M. Halberstadt, F. Beby, F. Flueckiger, M. Boehnke and J. G. Garweg. 2005. "Toxic effects of indocyanine green, infracyanine green, and trypan blue on the human retinal pigmented epithelium." *Graefes Archive for Clinical and Experimental Ophthalmology* 243(9): 917-925.
47. Veckeneer, M., K. Overdam, J. Monzer, K. Kobuch, W. Marle, H. Spekrijse and J. Meurs. 2001. "Ocular toxicity study of trypan blue injected into the vitreous cavity of rabbit eyes." *Graefes Archive for Clinical and Experimental Ophthalmology* 239(9): 698-704.
48. Gonzales-Tello, P., F. Camacho and G. Blazquez. 1994. "Density and Viscosity of Concentrated Aqueous Solutions of Polyethylene Glycol." *Journal of Chemical & Engineering Data* 39(3): 611-614.

Vioron®

The versatile trypan blue dye for the anterior segment



Brilliant visualization of the anterior lens capsule

Excellent distinction of the capsulorhexis margin

Quick and easy application (Ready-to-use)

Approved for DMEK

Composition and properties


1 ml of Vioron® contains:

- 0.6 mg Trypan Blue
- Disodium phosphate ($\text{Na}_2\text{HPO}_4 \times 2 \text{H}_2\text{O}$)
- Monosodium phosphate ($\text{NaH}_2\text{PO}_4 \times 2 \text{H}_2\text{O}$)
- Sodium chloride (NaCl)
- Water for injection

Density: 1.00 g/cm³


Packaging units

 **G-81002 Vioron® Syringe**
0.5 ml syringe, 5 pcs. per box, sterile

 **G-81001 Vioron® Vial**
0.5 ml vial, 5 pcs. per box, sterile

Fields of application

Vioron® was developed for ophthalmic surgical procedures in the anterior segment of the eye such as cataract operations or keratoplasties. Staining the anterior lens capsule makes it more visible, thus facilitating capsulorhexis and minimizing the risk of tearing. Furthermore, Vioron® facilitates the preparation and transfer of the donor cornea in the case of lamellar corneal transplantations and the removal of the diseased Descemet's membrane in case of DMEK and DS(A)EK).

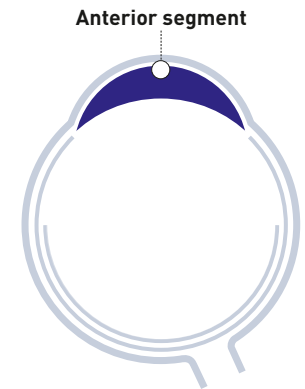

References

Keratoplasty

49. Melles, G.R.J., T. San Ong, B. Ververs and J. van der Wees. 2006. "Descemet Membrane Endothelial Keratoplasty (DMEK)." *Cornea* 25(8): 987-990. 50. Ham L., I. Dapena, C. van Luijk, J. van der Wees and G. R. J. Melles. 2009. "Descemet Membrane Endothelial Keratoplasty (DMEK) for Fuchs endothelial dystrophy: review of the first 50 consecutive cases." *Eye* 23(10): 1990-1998. 51. Cursiefen, C. and F. E. Kruse. 2010. "DMEK: Posteriore lamelläre Keratoplastik-technik." *Der Ophthalmologe* 107(4): 370-376. 52. Droutsas, K., L. Ham, I. Dapena, G. Geerling, S. Oellerich and G.R.J. Melles. 2010. "Visus nach Descemet-Membran Endothelkeratoplastik (DMEK): Ergebnisse der ersten 100 Eingriffe bei Fuchs'scher Endotheldystrophie." *Klinische Monatsblätter für Augenheilkunde* 227(06): 467-477. 53. Dapena, I., K. Moutsouris, K. Droutsas, L. Ham, K. van Dijk and G.R.J. Melles. 2011. "Standardized "To-Touch" Technique for Descemet Membrane Endothelial Keratoplasty." *Archives of Ophthalmology* 129(1): 88-94. 54. Dirisamer, M., L. Ham, I. Dapena, K. Moutsouris, K. Droutsas, K. van Dijk, L. E. Frank, S. Oellerich and G.R.J. Melles. 2011. "Efficiency of Descemet Membrane Endothelial Keratoplasty: Clinical Outcome of 200 Consecutive Cases after a Learning Curve of 25 Cases." *Archives of Ophthalmology* 129(11): 1435-1443. 55. Heindl, L. M., S. Riss, B. O. Bachmann, K. Laaser, F. E. Kruse and C. Cursiefen. 2011. "Split cornea transplantation for 2 recipients. A new strategy to reduce corneal tissue cost and shortage." *Ophthalmology* 118(2): 294-301. 56. Güell, J. L., M. Morral, O. Gris, D. Elies and F. Manero. 2013. "Bimanual technique for insertion and positioning of Endothelial-Descemet Membrane graft in Descemet Membrane Endothelial Keratoplasty." *Cornea* 32(12): 1522-1526. 57. Ciechanowski, P. P., K. Droutsas, L. Baydoun, M. Dirisamer, S. Oellerich, G.R.J. Melles. 2014. "Standardisierte Descemet-Membran-Endothelkeratoplastik (DMEK). Technik und aktuelle Ergebnisse." *Der Ophthalmologe* 111(11): 1041-1049. 58. Baydoun, L., K. van Dijk, I. Dapena, F. U. Musa, V. Liarakos, L. Ham and G.R.J. Melles. 2015. "Repeat Descemet Membrane Endothelial Keratoplasty after complicated primary Descemet Membrane Endothelial Keratoplasty." *Ophthalmology* 122(1): 8-16.

Capsulorhexis



59. Melles, G.R.J., P. W. T. de Waard, J. H. Pameijer and W. H. Beekhuis. 1999. "Trypan blue capsule staining to visualize the capsulorhexis in cataract surgery." *Journal of Cataract and Refractive Surgery* 25(1):7-9. 60. Melles, G.R.J., P. W. T. de Waard, J. H. Pameijer and W. H. Beekhuis. 1999. "Färbung der Linsenkapsel mit Trypanblau zur Visualisierung der Kapsulorhexis bei Maturkataraktchirurgie." *Klinische Monatsblätter für Augenheilkunde* 215(12): 342-344. 61. Werner, L., S. K. Pandey, M. Escobar-Gomez, D. S. Hoddinott and D. J. Apple. 2000. "Dye-enhanced cataract surgery. Part 2: Learning critical steps of phacoemulsification1." *Journal of Cataract and Refractive Surgery* 26(7): 1060-1065. 62. De Waard, P.W.T., C. J. Budo and G.R.J. Melles. 2002. "Trypan blue capsular staining to "find" the leading edge of a "lost" capsulorhexis." *American Journal of Ophthalmology* 134(2): 271-272. 63. Jacob, S., A. Agarwal, A. Agarwal, S. Agarwal, S. Chowdhary, R. Chowdhary and A. A. Bagmar. 2002. "Trypan blue as an adjunct for safe phacoemulsification in eyes with white cataract." *Journal of Cataract and Refractive Surgery* 28(10): 1819-1825. 64. Özkiris, A., O. Arslan, E. Cıcık, N. Köylüoğlu and C. Evereklioglu. 2003. "Open-sky capsulorhexis in triple procedure: with or without trypan blue?" *European Journal of Ophthalmology* 13(9-10): 764-769. 65. Nardi, M., U. Benelli, M. Figus and M. P. Bartolomet. 2004. "Open-sky capsulorhexis in triple procedure: with or without trypan blue?" *European Journal of Ophthalmology* 14(4): 358-358. 66. Singh, A. J., U. A. Sarodia, L. Brown, R. Jagjivan and R. Sampath. 2003. "A histological analysis of lens capsules stained with trypan blue for capsulorhexis in phacoemulsification cataract surgery." *Eye* 17(5): 567-570. 67. Xiao, Y., Y. H. Wang, Z. Y. Fu and H. Hong. 2004. "Staining the anterior capsule with indocyanine green or trypan blue for capsulorhexis in eyes with white cataract." *International Ophthalmology* 25(5-6): 273-276. 68. Chung, C. F., C. C. Liang, J. S. Lai, E. S. Lo and D. S. Lam. 2005. "Safety of trypan blue 1% and indocyanine green 0.5% in assisting visualization of anterior capsule during phacoemulsification in mature cataract." *Journal of Cataract and Refractive Surgery* 31(5): 938-942. 69. Wong, V. W. Y., T. Y. Y. Lai, G. K. Y. Lee, P. T. H. Lam and D. S. C. Lam. 2006. "A prospective study on trypan blue capsule staining under air vs under viscoelastic." *Eye* 20(7): 820-825. 70. Chéour, M., F. B. Brahim, A. Zarrad, N. Khémiri, K. Mghaieth and A. Kraiem. 2007. "Phacoémulsification des cataractes blanches en utilisant le bleu trypan." *Journal Français d'Ophtalmologie* 30(9): 914-917. 71. Ziakas, N. G., K. Boboridis, E. Nakos, D. Mikropoulos, V. Margaritis and A.G.P. Konstas. 2009. "Does the use of trypan blue during phacoemulsification affect the intraocular pressure?" *Canadian Journal of Ophthalmology* 44(3): 293-296. 72. Smith, E. F., R. U. Desai, A. Schrier, B. Enriquez and B. K. Purewal. 2010. "Trypan blue capsulorhexis." *Ophthalmology* 117(7): 1462-1462.e1. 73. Jaber, R., L. Werner, S. Fuller, S. C. Kavoussi, S. McIntyre, M. Burrow and N. Mamalis. 2012. "Comparison of capsulorhexis resistance to tearing with and without trypan blue dye using a mechanized tensile strength model." *Journal of Cataract and Refractive Surgery* 38(3): 507-512. 74. Giammaria, D., M. Giannotti, A. Scopelliti, G. Pellegrini and B. Giannotti. 2013. "Under-air staining of the anterior capsule using trypan blue with a 30 g needle." *Clinical Ophthalmology [Auckland, NZ]*, 7: 233-235.

F-Octane F-Decalin

Ultrapure perfluorocarbons
for intraoperative tamponades




- Gentle retinal unfolding and stabilization
- Drainage of subretinal fluids
- Refloating luxated lenses
- Short-term tamponade
- Outstanding stability and biocompatibility
- Ready-to-use syringes

Packaging units		F-Octane	F-Decalin
	Syringe 5 ml	G-80315	G-80115
	Syringe 7 ml	G-80317	G-80117
	Vial 5 ml	G-80305	G-80105
	Vial 7 ml	G-80307	G-80107

Field of application

F-Octane and F-Decalin are used as medical adjuvants for gentle retinal unfoldings, giant tears, traumata, laser coagulation as well as cryotherapy. Furthermore, they are used for refloating luxated lenses and as short-term tamponades.

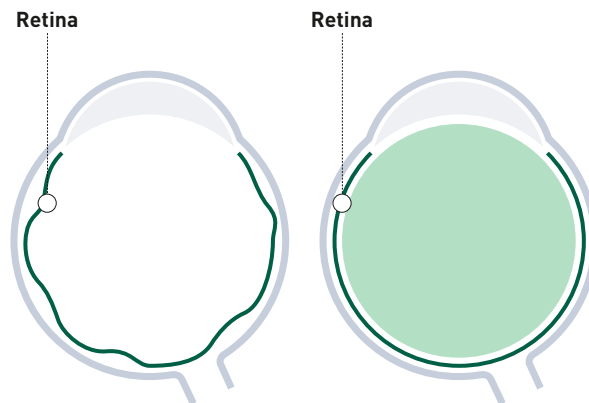

Composition and properties

F-Octane and F-Decalin are sterile fluorocarbon compounds with high density (1.76 g/cm³ and 1.93 g/cm³). They only consist of C-C and C-F bonds and do not contain any relevant amounts of biologically active components, thanks to the complex purification process at Fluoron. Due to the exceptional strength of the C-F bonds, F-Octane and F-Decalin are chemically and physiologically inert and absolutely non-toxic.

	F-Octane	F-Decalin
Density [g/cm³]	1.76	1.93
Vapor pressure [mbar] at 25° C	18.5	8.0
Refractive Index at 20° C	1.2700	1.3110
Surface tension [mN/m] at 25° C	14.0	19.0
Interface tension [mN/m] at 25° C	55.0	57.8
Composition	100 % Perfluorooctan (PFO)	100 % Perfluorodecalin (PFD)

References

75. Chang, S., V. Reppucci, N. J. Zimmerman, M. H. Heinemann and D. J. Coleman. 1989. "Perfluorocarbon liquids in the management of traumatic retinal detachments." *Ophthalmology* 96(6): 785-792. 76. Meinert, H. 1992. "Requirements of perfluorocarbons for use in ophthalmology." *Journal of Vitreoretina* 1, 5-16. 77. Binder, S., M. Velikay, A. Wedrich, U. Stolba and P. Datlinger. 1992. "Die klinische Anwendung flüssiger Perfluorocarbone in der Netzhautchirurgie." *Spektrum Augenheilkunde* 6(1): 4-7. 78. Le Mer, Y. and P. Kroll. 1991. "Flüssiges Perfluorocarbon bei der Behandlung von Riesenrissen." *Klinische Monatsblätter für Augenheilkunde* 199(10): 256-258. 79. Chang, S., H. Lincoff, N. J. Zimmerman and W. Fuchs. 1989. "Giant Retinal Tears Surgical techniques and results using perfluorocarbon liquids." *Archives of Ophthalmology* 107(5): 761-766.

EasyGas[®]

The first ready-to-use gas tamponade



Quick and easy application through sterile, pre-filled system

Safe usage because of precise, non-expanding mixture ratio in each syringe

No mix-up of gases due to colour coding

Three gases for different tamponade durations


Reduced risk for hypertension or ischemia, because manual mixing is not required


No subsequent surgery for tamponade removal necessary

Contains patient information card and patient wristband

Packaging units


 **G-80950 EasyGas[®] SF₆**
Syringe 40 ml, sterile

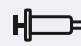 **G-80960 EasyGas[®] C₂F₆**
Syringe 40 ml, sterile

 **G-80970 EasyGas[®] C₃F₈**
Syringe 40 ml, sterile

Fields of application

EasyGas® SF₆, EasyGas® C₂F₆ and EasyGas® C₃F₈ are the first ready-to-use gas tamponades. The sterile, pre-filled, ready-to-use system offers a quick and easy application of the tamponades. EasyGas® is used as long-term tamponade after operative treatment of severe retinal detachment.

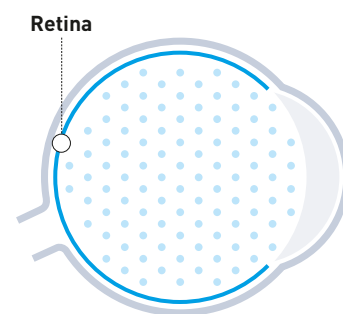

	EasyGas® SF ₆	EasyGas® C ₂ F ₆	EasyGas® C ₃ F ₈
Effective tamponade time [days]	6	15	30
Retention time / Longevity* [weeks]	1–2	4–5	6–8
Non-expansive gas concentration* [%]	20	16	12

Indication

	EasyGas® SF ₆	EasyGas® C ₂ F ₆	EasyGas® C ₃ F ₈
Retinal detachment with giant tears	✓	✓	–
Retinal detachment without proliferation	✓	✓	–
Retinal detachment in cases of proliferative diabetic retinopathy (PDR)	✓	✓	–
Proliferative vitreoretinopathy (PVR)	✓	✓	✓
Traumatic retinal detachment	✓	✓	–

Reference: Cekic O., Ohji M. (2010): Intraocular Gas Tamponades. Semin Ophthalmol. 2000 Mar; 15 (1): 3-14.

Composition and properties

	EasyGas® SF ₆	EasyGas® C ₂ F ₆	EasyGas® C ₃ F ₈
Density of pure gas	6.17 kg/m ³	5.73 kg/m ³	8.17 kg/m ³
Density of ready-to-use mixture	2.34 kg/m ³	1.95 kg/m ³	2.01 kg/m ³
Composition	20 % SF ₆ 5.0 80 % synthetic air 6.0	16 % C ₂ F ₆ 5.0 84 % synthetic air 6.0	12 % C ₃ F ₈ 4.0 88 % synthetic air 6.0

Siluron®

A new generation of silicone oils
with an innovative molecular design



For incision
sizes 23 + 25
Gauge

The new generation of Siluron® silicone oils is characterized by its special property of a significantly higher emulsification resistance. This is based on an intelligent mixture of different chain lengths of molecules and the resulting dynamic viscosity. Good injectability in cases of small incisions is a further advantage of these innovative silicone oils.

High resistance to emulsification

Short injection time

Good long-term tolerability

Excellent chemical purity


Siluron® 2000


Premium silicone oil with
customized extensional viscosity

Siluron® XTRA

Premium silicone oil with
an Xtra portion of elasticity

Packaging units

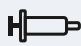 **G-80740 Siluron® 2000 Syringe**
10 ml, sterile


 **G-80750 Siluron® Xtra Syringe**
10 ml, sterile

Fields of application

Siluron® silicone oils are used as long-term tamponades after operative treatment of severe retinal detachment, particularly for:

- Retinal detachments with giant tears
- Retinal detachments with proliferative vitreoretinopathy (PVR)
- Retinal detachments in cases of proliferative diabetic retinopathy (PDR)
- Traumatic retinal detachments

Due to their specific density of 0.97g/cm³ the Siluron® silicone oils float on water.

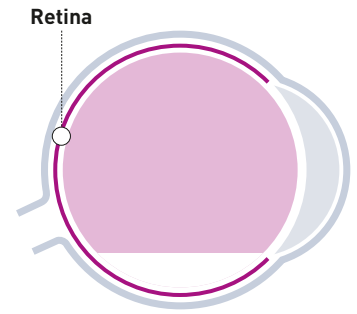

Good long-term tolerability

Excellent chemical purity

Chemically and physiologically inert

Siluron® 1000


Easily injectable

Siluron® 5000

Higher resistance to emulsification

Packaging units

 **G-80720 Siluron® 1000 Syringe** 10 ml, sterile
G-80820 Siluron® 5000 Syringe 10 ml, sterile

 **G-80710 Siluron® 1000 Vial** 10 ml, sterile
G-80810 Siluron® 5000 Vial 10 ml, sterile

References

80. Neuer, K. L., S. Bohnacker, N. Feucht, C. P. Lohmann and M. Maier. 2017. "In Vitro Biocompatibility of Silicone Oil Siluron Xtra on Porcine Retina in a Perfusion Culture System." *Journal of Clinical and Experimental Ophthalmology* 8(679):2. 81. Hussain, R. N., J. Myneni, T. Stappler and D. Wong. 2017. "Polydimethyl Siloxane as an Internal Tamponade for Vitreoretinal Surgery." *Ophthalmologica* 238(1-2): 68-73. 82. Shen, Y.-D. and C.-M. Yang C.-M. 2007. "Extended silicone oil tamponade in primary vitrectomy for complex retinal detachment in proliferative diabetic retinopathy: A long-term follow-up study" *European Journal of Ophthalmology* 17 (6): 954-960. 83. Lucke, K. H., M. H. Foerster and H. Laqua. 1987. "Long-Term results of vitrectomy and silicone oil in 500 cases of complicated retinal detachments." *American Journal of Ophthalmology* 104(6): 624-633. 84. Heidenkummer, H. P., A. Kampik and S. Thierfelder. 1991 "Emulsification of silicone oils with specific physicochemical characteristics." *Graefe's Archive for Clinical and Experimental Ophthalmology* 229(1): 88-94. 85. Petersen, J. 1987. "The physical and surgical aspects of silicone oil in the vitreous cavity." *Graefe's Archive for Clinical and Experimental Ophthalmology* 225(6): 452-456. 86. Han, L., J. D. Cairns, W. G. Campbell, M. F. McCombe, W. J. Heriot and J. B. Heinze. 1998. "Use of silicone oil in the treatment of complicated retinal detachment: results from 1981 to 1994." *Australian and New Zealand Journal of Ophthalmology*

26(4), 299-304. 87. Kampik, A., V. P. Gabel and D. Spiegel. 1984. "Intraokulare Tamponade mit hochviskosem Silikonöl bei massiver proliferativer Vitreo-Retinopathie." *Klinische Monatsblätter für Augenheilkunde* 185(11): 368-370. 88. Stalmans, P., A. M. Pinxten and D. S. Wong. 2015. "Cohort safety and efficacy study of siluron2000 emulsification-resistant silicone oil and F4H5 in the treatment of full-thickness macular hole." *Retina* 35(12): 2558-2566. 89. Caramoy, A., N. Hagedorn, S. Fauser, W. Kugler, T. Gross and B. Kirchhof. 2011. "Development of emulsification-resistant silicone oils: can we go beyond 2000 mPas silicone oil? *Investigative ophthalmology & visual science* 52(8): 5432-5436. 90. Caramoy, A., S. Schröder, S. Fauser and B. Kirchhof. 2010. "In vitro emulsification assessment of new silicone oils." *British Journal of Ophthalmology* 94(4): 509-512. 91. Chan, Y. K., C. O. Ng, P. C. Knox, M. J. Garvey, R. L. Williams, and D. Wong. 2011. "Emulsification of silicone oil and eye movements." *Investigative ophthalmology & visual science* 52(13): 9721-9727. 92. Williams, R. L., M. J. Day, M. J. Garvey, G. Morphis, C. Irigoyen, D. Wong and T. Stappler. 2011. "Injectability of silicone oil-based tamponade agents." *British Journal of Ophthalmology* 95(2): 273-276. 93. Williams, R. L., M. Day, M. J. Garvey, R. English and D. Wong, D. 2010. "Increasing the extensional viscosity of silicone oil reduces the tendency for emulsification." *Retina* 30(2): 300-304.


Overview of properties

Physicochemical properties of Siluron® oils

Property	Siluron® 1000	Siluron® 5000	Siluron® 2000	Siluron® XTRA
Density [g/cm ³] at 25° C	0.97	0.97	0.97	0.97
Viscosity [mPas] at 25° C	900 - 1200	4800 - 5500	2000 - 2400	4100 - 4800
Refractive index	1.404	1.404	1.404	1.404
Solubility in water	non miscible	non miscible	non miscible	non miscible
Composition [w%]	100 % Polydi-methylsiloxane (PDMS)	100 % Polydi-methylsiloxane (PDMS)	95 % Siluron® 1000 + 5 % PDMS (2.5 M mPas)	90 % Siluron® 1000 + 10 % PDMS (2.5 M mPas)
Elasticity (Je ⁰) [Pas]	2 x 10 ⁻⁵	1 x 10 ⁻⁵	6.5 x 10 ⁻⁴	1.4 x 10 ⁻³
Shear viscosity (at 8.37 s ⁻¹ and 37° C) [mPas]	931	4303	1800	4377
Volatile components (200° C, 24 h) [%]	≤ 0.2 %	≤ 0.2 %	≤ 0.2 %	≤ 0.2 %

Reference: Caramoy A., Hagedorn N., Fauser S., Kugler W., Gross T., Kirchhof B.: Development of emulsification-resistant silicone oils: can we go beyond 2000 mPas silicone oil? Invest Ophthalmol Vis Sci. 2011; 52: 5432-5436

Comparison of emulsification rate

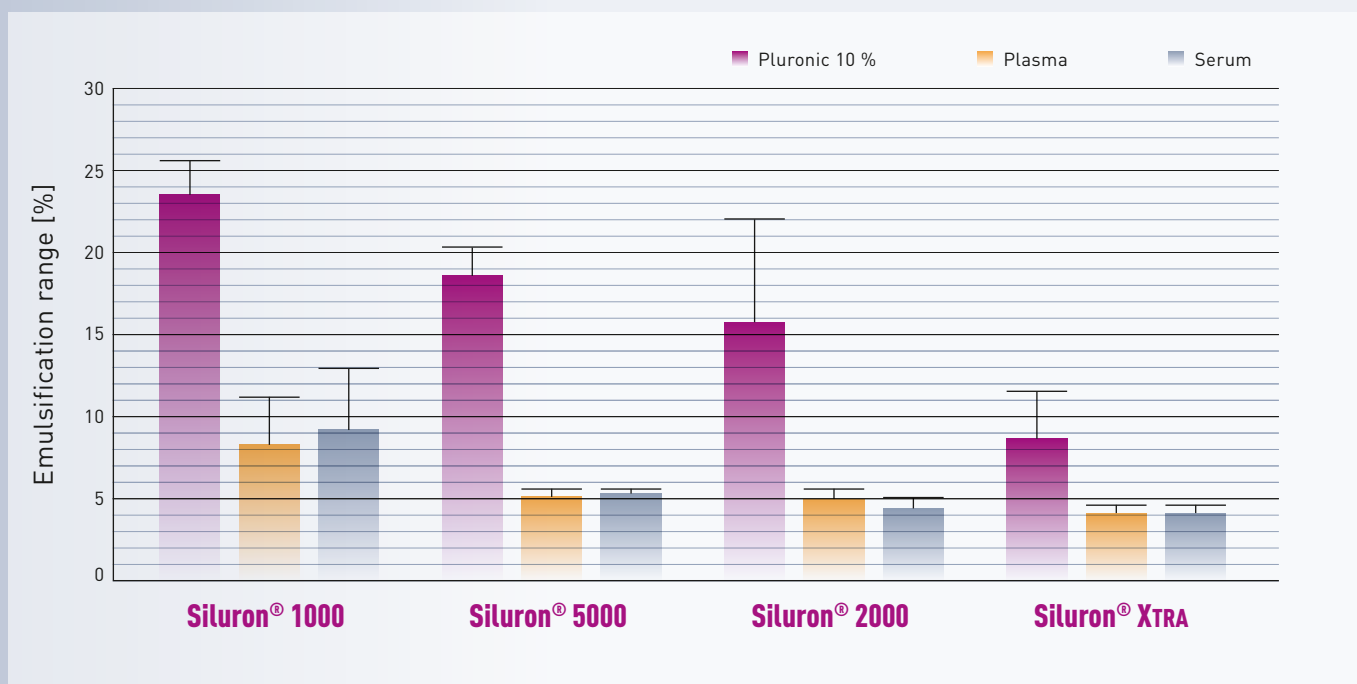


Reference: Caramoy A., Hagedorn N., Fauser S., Kugler W., Gross T., Kirchhof B.: Development of emulsification-resistant silicone oils: can we go beyond 2000 mPas silicone oil? Invest Ophthalmol Vis Sci 2011; 52: 5432-5436

Comparison of emulsification resistance



Reference: Wong et. al

Comparison of injection time 5,5 bar injection pressure, 20 gauge injection cannula

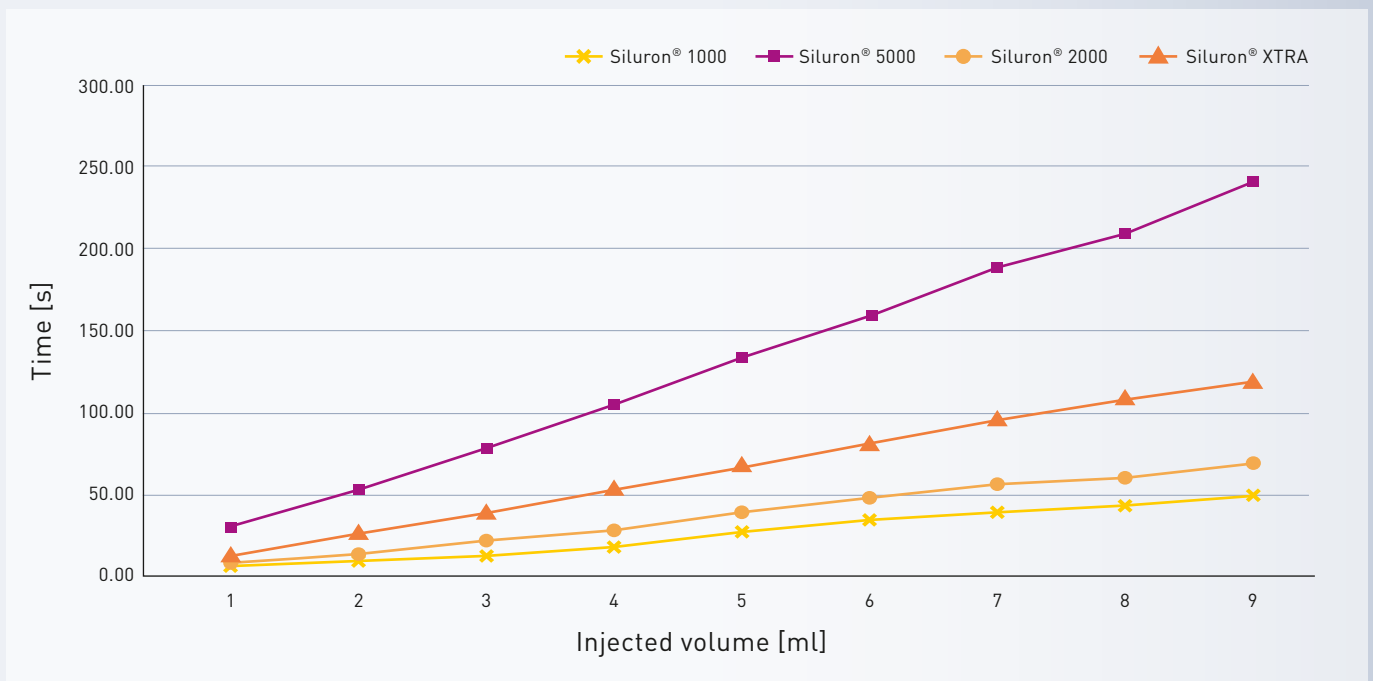

Reference: Williams RL., Day MJ., Garvey MJ., Morphis G., Irigoyen C., Wong D., Stappler T.: Injectability of silicone oil-based tamponade agents. Br J Ophthalmol. 2011; 95: 273-276

Densiron® XTRA

The heavy silicone oil
with molecular design



Unique molecularly design

Heavier than water

Ideal for inferior pathologies

Removing proliferative milieu
in lower part of retina

Avoiding unpleasant constraints
for patient ("head-down-position")

Easy to inject

25G compatible

Highly resistant against emulsification

Caramoy A., Schröder S., Fauser S., Kirchhof B. (2010) In vitro emulsification assessment of new silicone oils. Br J Ophthalmol 94, 509-512

Composition and properties

Density [g/cm³] at 25° C:

1.05 - 1.07

Viscosity [mPas] at 25° C:

1.000 - 1.400

Composition [w%]:

30.5 % F6H8

69.5 % Siluron® Xtra

Packaging units




G-80925 Densiron® Xtra Syringe

10 ml syringe, 1 pc. per box, sterile


Fields of application

Densiron® Xtra is used as a temporary tamponade after operative treatment of severe retinal detachment, particularly for:

- Inferior and posterior retinal holes
- Retinal detachments with giant tears
- Retinal detachments with proliferative vitreoretinopathy (PVR)
- Retinal detachments in cases of proliferative diabetic retinopathy (PDR)
- Traumatic retinal detachments

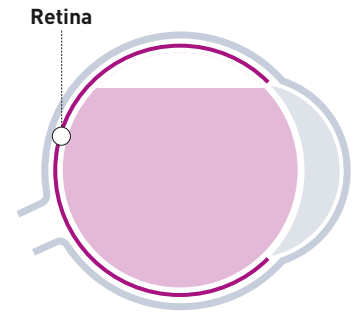


Injection times Densiron® Xtra

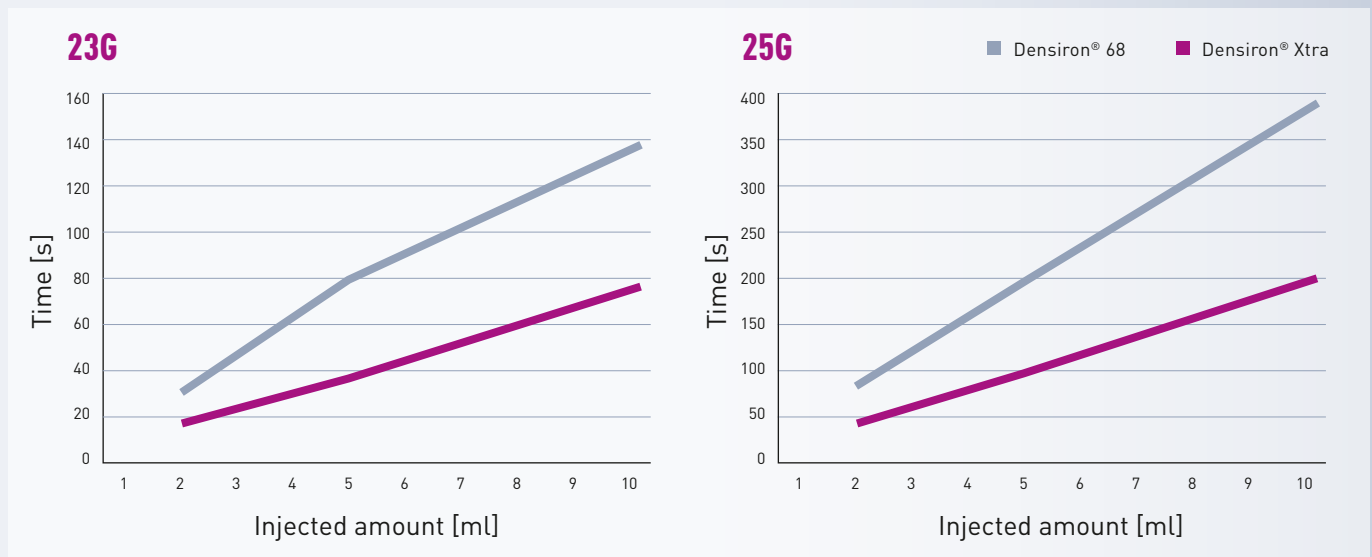

Vitrectomy System, 6 bar

In vitro emulsification ranges of various silicone oils when using plasma as emulsifier

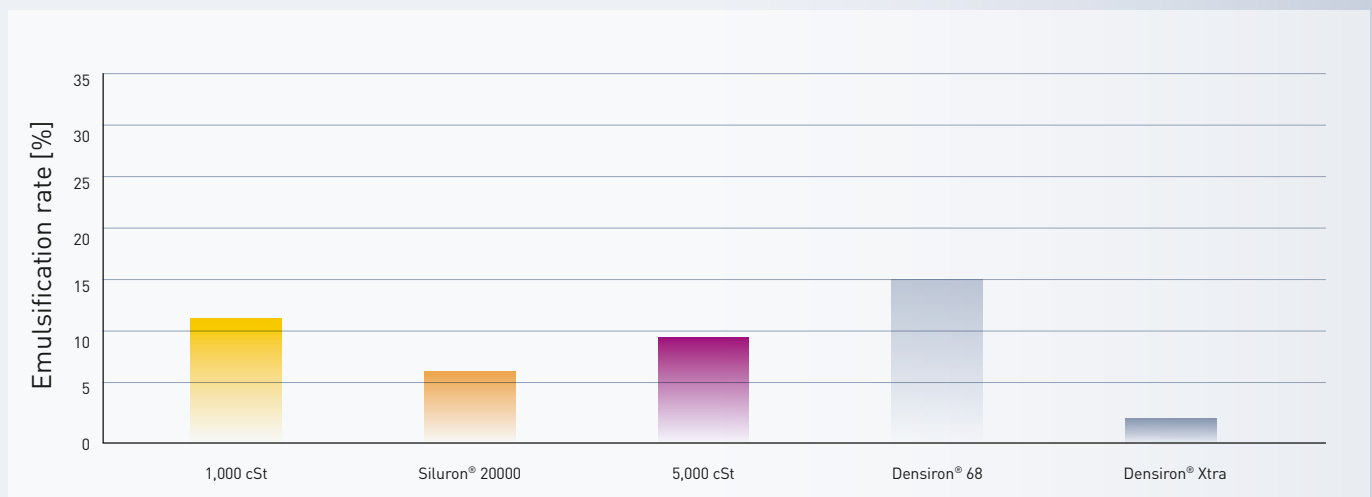

Reference: "In vitro emulsification assessment of new silicone oils." Caramoy et al. Br J Ophthalmol. 2010 Apr;94(4):509-12

Testimonials Densiron® Xtra

“Temporal inverted ILM flap technique combined with heavy silicone oil (Densiron Xtra) for macular detachment associated with ODP is a highly effective alternative technique. This procedure achieved very rapid resolution of the submacular fluid with successful anatomical and functional results.”

Oncel, M: A Novel Approach for the Management of Macular Detachment Associated with Optic Disc Pit: Temporal Inverted Internal Limiting Membrane Flap Technique and a New Heavy Silicone Oil (Densiron Xtra)



Prof. Francesco Boscia

MD, Associate Professor and Chair at the Department of Ophthalmology at the Sassari University, Sassari, Sardegna (IT)

Which are the key pathologies and why?

“As tamponade in recurrent inferior rhegmatogenous retinal detachment, especially if complicated by severe proliferative vitreoretinopathy.”

What features do you like most?

“It can effectively tamponade inferior retina with the patient standing upright (...). I routinely use 25G system and I never met any trouble in injecting and aspirating Densiron Xtra.”

What is your conclusion about Densiron Xtra?

“It’s an essential surgical tool for every vitreoretinal surgeon who needs to face with complex pathologies. It is effective in tamponing and stabilising the inferior retina and safe at the same time.”



Dr. Vignesh Raja

Joondalup Eye Clinic and Perth Eye Hospital
Perth, Australia

Which are the key pathologies and why?

“I prefer to use Densiron Xtra for pathologies such as persistent macular hole, inferior retinal detachment with PVR, inability to posture face down, recurrent and chronic retinal detachment that need long term silicone oil endotamponade.”

What features do you like most?

“I like Densiron Xtra because of its heavier than water property, low risk of emulsification and low risk of developing retinal/macular toxicity. Removal of Densiron is straight forward (with the correct technique) with low risk of residual silicone oil bubbles.”

What is your conclusion about Densiron Xtra?

“Densiron Xtra adds to my retinal armamentarium and is my preferred agent for endotamponade in challenging and complicated cases.”



Dr. Theodor Stappler
Médecin adjoint, Unité de chirurgie
vitréorétinienne, Hôpital ophtalmique
Jules-Gonin, Lausanne (CH)

Which are the key pathologies and why?

"The treatment of inferior proliferation in recurrent retinal detachment (...) to exclude the aqueous environment containing cytokines and proliferative agents entirely from the retinal area which had just been treated."

What features do you like most?

"I can use Densiron Xtra irrespective of the gauge. The process of injection and removal has stopped being lengthy and arduous."

What is your conclusion about Densiron Xtra?

"Easy to inject and aspirate, decreased emulsification rate, yet heavy tamponade agent."



Dr. Andreas Kölbl
Ophthalmic Specialist,
Ophthalmic Private Practice, Eggenburg (AT)

Which are the key pathologies and why?

"Mainly for complicated retinal detachments (PVR) with tear formation and tensions in the inferior segment, also for tractions due to diabetic retinopathy and I'm happy with the results."

What features do you like most?

"The comparable low viscosity and hence the excellent injectability even via 25G systems (...)."

What is your conclusion about Densiron Xtra?

"I use Densiron Xtra because I feel more secure in complicated retinal detachments with pathologies in the inferior segment for elderly patients for whom correct patient positioning cannot be guaranteed."



**Dr. Antonio
Palomino Muñoz**
Oftalmólogo, Hospital Quiron San Jose,
Madrid (ES)

Which are the key pathologies and why?

"We use it in all retinal surgery in which are predisposing factors for PVR."

What features do you like most?


"The quality that I appreciate most is the ease of injection, even with 25G. Also its intraocular tolerance and stability against emulsification is appreciable."

What is your conclusion about Densiron Xtra?

"These referred qualities make Densiron Xtra an important ally in complex vitreoretinal surgery improving its prognosis."

F4H5[®] WashOut

The simple solution for oil residues in vitreoretinal surgery

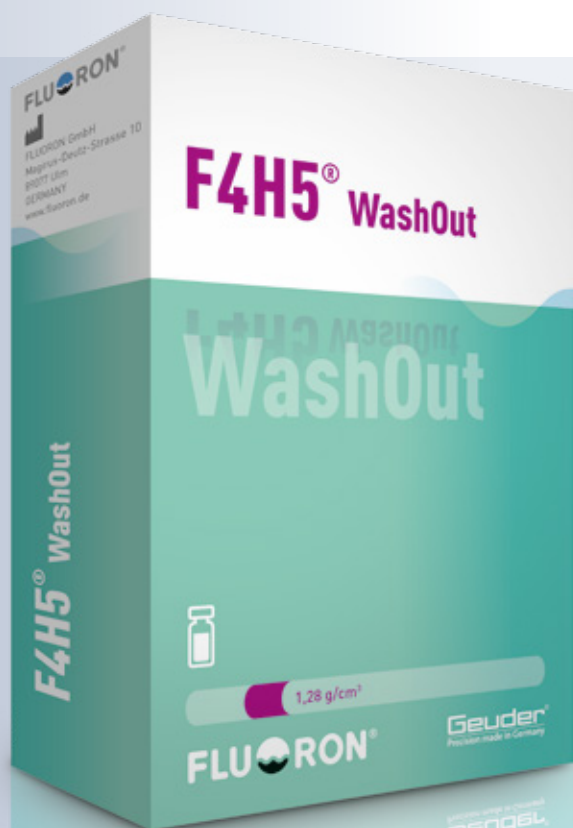

Unique amphiphilic properties

Solves silicone oil efficiently

Removes silicone oil residues and "sticky oil"

Rinses silicone oil-polluted IOL

Biocompatible

Composition and properties

Density [g/cm³] at 25° C: 1.28

Viscosity [mPas] at 25° C: 1.05

Mix ratio F4H5 : Silicone oil:

Mix in all ratios

Packaging units



G-80615 F4H5[®] WashOut Vial

5 ml vial, 1 pc. per box, sterile

Fields of application

F4H5[®] WashOut is a biocompatible solvent for removing silicone oil residues from the retina and for cleaning intraocular lenses after silicone oil tamponades.

Video

Scan QR-Code for further information on F4H5[®] WashOut



References

94. Stappler, T., R. Williams and D. Wong. 2010. "F4H5 - A novel substance for the removal of silicone oil from intraocular lenses." *British Journal of Ophthalmology* 94(3): 364-367.
95. Chan, Y. K., H. C. Cheng, J. Wu, Y. H. M. Tang, S. T. Chan, D. Wong and H. C. Shum, H. 2018. "A perfluorobutylpentane (F4H5)-based solution for the removal of residual emulsified silicone oil." *Acta ophthalmologica* 96(1): e38-e45.
96. Coppola, M., C. Del Turco, G. Querques and F. Bandello. 2017. "Perfluorobutylpentane (F4H5) solvent-assisted silicone oil removal technique." *Retina* 37(4): 793-795.
97. Stalmans, P. and D. Wong. 2015. "Cohort Safety & efficacy

study of Siluron 2000 emulsification-resistant silicone oil and F4H5 in the treatment of full-thickness macular hole." *Retina* 35(12): 2558-2566.
98. Tatsios, J., W. Kugler, A. M. Jousseaume. 2016. "Problem "gelöst": F4H5 als Auswaschlösung für Silikonreste in der Glaskörperchirurgie." *German Medical Science GMS Publishing House Doc16rg28*.
99. Liang, Y., N. Kociok, M. Leszczuk, W. Hiebl, B. Theisinger, A. Lux and A. M. Jousseaume. 2008. "A cleaning solution for silicone IOLs and 'sticky silicone oil'." *British Journal of Ophthalmology* 92: 1-1.

F6H8[®] Vitreous Substitute

The rinsing tamponade
in vitreoretinal surgery



Temporary endotamponade
in cases of complicated, especially
inferior retinal detachments

“Third hand“ when relocating
the macula

Retinal unfolding with
low contact pressure

Rinses IOLs after
silicone oil tamponades

Composition and properties

Density [g/cm³] at 25° C: 1.33

Viscosity [mPas] at 25° C: 3.44

Mix ratio F6H8 : Silicone oil:

optimal 50:50 to 30:70

Packaging units



G-80606 Vitreous Substitute F6H8[®] Vial
6 ml vial, 1pc. per box, sterile



G-80609 Vitreous Substitute F6H8[®] Vial
9 ml vial, 1 pc. per box, sterile

References

100. Wong, D and N. Lois. 2000. "Perfluorocarbons and Semifluorinated Alkanes." Seminars in Ophthalmology Vol. 15, No. 1: 25-35. Taylor & Francis. 101. Meinert, H. and T. Roy. 2000. "Semifluorinated alkanes – A new class of compounds with outstanding properties for use in ophthalmology." European Journal of Ophthalmology 10(3): 189-197. 102. Schatz, B., Y. El-Shabrawi, A. Haas and G. Langmann. 2004. "Adverse side effects with Perfluorohexyloctane as a long-term tamponade agent in complicated vitreoretinal surgery." Retina 24(4): 567-573.

Fields of application

F6H8[®] serves as an endotamponade in complicated retinal detachments and as an intraoperative tool for retinal surgery. Due to its significantly lower density compared to conventional perfluorocarbon liquids, F6H8[®] offers enormous advantages particularly for retinal translocation. F6H8[®] is also an excellent biocompatible solvent for the removal of silicone oil residues from the vitreous chamber as well as being suitable for cleaning intraocular lenses after a silicone oil tamponade.

103. Langefeld, S., B. Kirchof, H. Meinert, T. Roy, A. Aretz and N. F. Schrage. 1999. "A new way of removing silicone oil from the surface of silicone intraocular lenses." Graefes' archive for clinical and experimental ophthalmology 237(3): 201-206. 104. Kobuch, K., D. H. Menz, H. Hoerauf, J. H. Dresch and V. P. Gabel. 2001. "New substances for intraocular tamponades: perfluorocarbon liquids, hydrofluorocarbon liquids and hydrofluorocarbon-oligomers in vitreoretinal surgery." Graefes' archive for clinical and experimental ophthalmology 239(9): 635-642.

Accessories*

Accessories PFCL




G-33057

Chang

PFCL Cannula


for injection of heavy fluids
dual bore, coaxial
25 gauge / 0.5 mm tip
20 gauge / 0.9 mm shaft

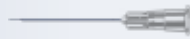

G-34285

Single-Use PFCL Cannula

for injection of heavy fluids
dual bore, coaxial
23 gauge / 0.64 mm
10 pcs. per box, sterile

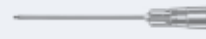

Single-Use

Backflush Cannula

with silicone brush
5 pcs. per box, sterile

G-34293 20 gauge / 0.9 mm

G-34294 23 gauge / 0.6 mm

G-34297 25 gauge / 0.5 mm



G-37002

Backflush Handpiece

with silicone chamber and
Luer-Lock connector
overall length 115 mm



G-34289

Single-Use

Backflush Handpiece

with silicone chamber, Luer-Lock
10 pcs. per box, sterile



Single-Use

Backflush Cannula

blunt tip

5 pcs. per box, sterile

G-34291 20 gauge / 0.9 mm

G-34296 23 gauge / 0.6 mm

G-34299 25 gauge / 0.5 mm

Accessories EasyGas®



G-80975

Single-Use Injection Cannula

for EasyGas®

30 gauge / 0.3 x 12 mm
100 pcs. per box, sterile



G-34492

Kirchhof

Single-Use Injection Cannula

for gas / viscous fluids
5.0 mm beveled tip with 4 infusion
side ports
2 metal sleeves, Luer-Lock plastic
adapter and 40 cm silicone tube
20 gauge / 0.9 mm
10 pcs. per box, sterile

* Accessories are optional. Further products available on demand.

Accessories silicone oils



Heidelberg Model
Cannula
 for injection or aspiration
 of viscous fluids
 and Densiron® 68, bevel 30°
G-32699 19 gauge / 1.1 mm
G-32698 18 gauge / 1.2 mm



G-33056
Roider
Aspiration Cannula
 for viscous fluids
 0.7 mm side port
 19 gauge / 1.0 mm



Single-Use Cannula
 to inject silicone oil
 10 pcs. per box, sterile
G-34497 20 gauge / 0.9 mm x 8 mm
G-34498 23 gauge / 0.6 mm x 8 mm



Hamburg Model
Injection Cannula
 for viscous fluids
 25 cm silicone tube with metal sleeve
 and Luer-Lock adapter
 3 spare silicone tubes
G-33470 20 gauge / 0.9 mm, beveled, 4 mm
G-33471 20 gauge / 0.9 mm, beveled, 5 mm
G-33472 20 gauge / 0.9 mm, beveled, 6 mm
G-33473 23 gauge / 0.6 mm, beveled, 4 mm
G-33474 23 gauge / 0.6 mm, beveled, 6 mm



Single-Use Injection Cannula
 for viscous fluids
 with 1 metal sleeve, Luer-Lock
 plastic adapter
 and 25 cm PVC tube
 20 gauge / 0.9 mm
 5 pcs. per box, sterile
G-33488 beveled tip, 4.0 mm
G-33489 beveled tip, 6.0 mm



G-33482
Single-Use Injection Cannula
 for viscous fluids, self-retaining,
 4.0 mm tip
 with 25 cm PVC tube
 Luer-Lock plastic adapter
 20 gauge / 0.9 mm
 5 pcs. per box, sterile



G-34480
Single-Use Injection Cannula
 for viscous fluids
 4.0 mm beveled tip
 with 2 metal sleeves, Luer-Lock
 plastic adapter
 and 25 cm silicone tube
 19 gauge / 1.1 mm
 10 pcs. per box, sterile

Accessories*

Accessories silicone oils



G-32697

Pressure Tube (reusable)
for injection of viscous fluids
Luer-Lock female / male



G-32696

Single-Use Pressure Tube
for injection of viscous fluid,
Luer-Lock female / male
10 pcs. per box, sterile



Stopper

for viscous fluid aspiration
with tube connection for single-use
syringe

G-33065 10 ml syringe

G-33066 20 ml syringe




G-28766

Single-Use Oil Injection System

to inject silicone oil pneumatically,
with protective cover for glass syringe,
pressure tube fits megaTRON® S3 / S4 HPS and
Pentasy®**, sterile

(G-28767 for megaTRON® and Accurus®, G-28768 for Millennium®)**

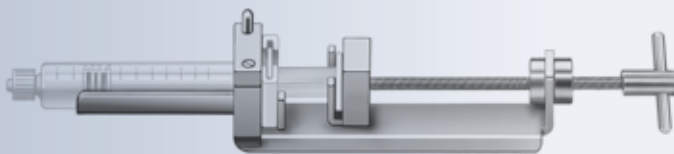

G-28752

Syringe Holder

for manual injection of
viscous fluids in glass syringes,
with clamp and retraction mechanism

* Accessories are optional. Further products available on demand.

** „megaTRON“ is a registered trademark of Geuder AG

„Pentasy“ is a registered trademark of Fritz RUCK Ophthalmologische Systeme GmbH

„ACCURUS“ is a registered trademark of Alcon Laboratories, Inc.

„Millennium“ is a registered trademark of MBI Millennium Biomedical, Inc.

Fluoron

Purity and variety made in Germany



Fluoron GmbH, based in Ulm, Germany, was founded in 1996 by Prof. Dr. Hasso Meinert and is a sister company of Geuder AG, Heidelberg. With his intellectual property rights, Prof. Meinert laid the foundation for a successful development of the company and accompanied Fluoron GmbH over 10 years on scientific topics. The company is managed by Mr Volker Geuder.

Fluoron GmbH develops and manufactures ultrapure innovative biomaterials for retinal and cataract surgery.

In this field, Fluoron GmbH plays a worldwide leading role in providing ophthalmic surgeons with creative and efficient solutions and consolidated its international competitive position by acquiring extensive intellectual property rights. The company's competence focuses on the development, manufacture and regulatory approval of light and heavy tamponades for retinal surgery, perfluorohydrocarbons and semifluorinated alkanes as temporary tamponades, as well as dyes for anterior and posterior segment surgery.

Trademarks:

„Brilliant Peel“, „Vioron“, „EasyGas“, „Siluron“, „Densiron“, „F4H5“ and „F6H8“ are registered trademarks of Fluoron GmbH.



Distributed in NZ by:

Medix 21

29 Jarden Mile Ngauranga 6035

Phone (04) 233 1627

Email admin@medix21.co.nz

Web www.medix21.co.nz

Fluoron GmbH reserves the right to make changes to compositions of their substances in response to recent developments. Fluoron GmbH and GEUDER AG do not assume liability for the accuracy of given statements. For questions regarding their products and availability please contact GEUDER AG.

Illustrations not drawn to scale.

Manufacturer: Fluoron GmbH Magirus-Deutz-Strasse 10 89077 Ulm Germany
Phone: +49 731 205 5997 0 Fax: +49 731 205 5997 28 info@fluoron.de www.fluoron.de

Distributor: GEUDER AG Hertzstrasse 4 69126 Heidelberg Germany
Phone: +49 6221 3066 Fax: +49 6221 303122 info@geuder.de www.geuder.de

Geuder[®]
Precision made in Germany

FLUORON[®]
PURITY AND VARIETY