

Distributed by:

Medix21 Surgical

29 Jarden Mile Ngauranga 6035

p 0800 633 4921 / 021 734 762

e camille@medix21.co.nz

w medix21surgical.co.nz



PRODUCT LISTING 2020
































Medicel AG
Dornierstrasse 11
9423 Altenrhein
Switzerland

Tel. +41 71 727 10 50
Fax +41 71 727 10 55
info@medicel.com
www.medicel.com

medicel
SWISS TECHNOLOGY FOR SURGERY

CONTENT

LENS INJECTION SYSTEM & CARTRIDGE SETS				PAGE
	ERGOJECT™	sterile	single-use top-loaded	 3
	ERGOJECT™	sterile	single-use back-loaded	 4
<i>NEW</i>	ACCUJECT™ DUAL	sterile	single-use top-loaded	 5
<i>NEW</i>	ACCUJECT™ DUAL	sterile	single-use back-loaded	 5
	ACCUJECT™	sterile	single-use top-loaded	 6
	ACCUJECT™	sterile	single-use back-loaded	 7
	ACCUJECT™ SCREW	sterile	single-use top-loaded	 8
	ACCUJECT™ SCREW	sterile	single-use back-loaded	 8
	VISCOJECT™ - BIO	sterile	single-use	 9
	VISCOJECT™ - BIO	non-sterile	reusable	 10
	VISCOJECT™ easy	sterile	single-use	 11,12
	VISCOJECT™ easy	non-sterile	reusable	 12
	VISCOJECT™ eco	sterile	single-use	 13
	MULTIJECT™	sterile	single-use	 14
	MULTIJECT™	non-sterile	reusable	 14
I/A INSTRUMENTS				
	Coaxial „Ergonomic design“	sterile	single-use	 15
	Coaxial „Soft-tip“	sterile	single-use	 15
	Bimanual „Ergonomic design“	sterile	single-use	 15
	Coaxial „Classical design“	sterile	single-use	 16
	Bimanual „Classical design“	sterile	single-use / reusable	 17
	Bimanual „Classical design“	non-sterile	reusable	 18
PHACO ACCESSORIES				
	4K F-Sonic Flared Phaco Sets	sterile	single-use	 19
	4K Phaco Accessories Sets	sterile, non-sterile	single-use, reusable	 20
	5R Phaco Accessories Sets	sterile, non-sterile	single-use, reusable	 21
	XXK Phaco Accessories Sets	sterile, non-sterile	single-use, reusable	 22
	COOLFLOW Infusion Sleeves	sterile, non-sterile	single-use, reusable	 23
	COMPATIBLES(Cassettes)	sterile	single-use	 24
DMEK INJECTOR SETS				
	ENDOJECT™	sterile	single-use	 25
	GÜELL DMEK	sterile	single-use	 25


ERGOJECT™

SINGLE-USE INJECTOR SETS, TOP-LOADED, STERILE

ERGOJECT™ is the world's first fully disposable screw injector with single-handed operation. Inspired by the Swiss watch industry, its precisely working micro gear provides an unprecedented smooth and controlled injection of almost any one-piece lens. Cartridges and top-loaded loading chambers are identical to the ones used in the ACCUJECT™ system. This makes the use of ERGOJECT™ very convenient for current ACCUJECT™ users.

TOP-LOADED / FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: ERGOJECT™ Injection system
with integrated cartridge and top-loaded loading chamber
(Box of 10)

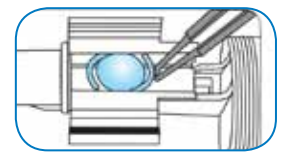

ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
ERT3000	ERGOJECT™ 3.0-TL Injector Set	2.4	2.7-2.8	3.1	1-piece lenses
ERT2600	ERGOJECT™ 2.6-TL Injector Set	2.2	2.5	2.8	1-piece lenses
ERT2400	ERGOJECT™ 2.4-TL Injector Set	2.1	2.4	2.7	1-piece lenses
ERT2200	ERGOJECT™ 2.2-TL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
ERT2100	ERGOJECT™ 2.1-TL Injector Set	1.9	2.1-2.2	2.4	1-piece lenses
ERT2000	ERGOJECT™ 2.0-TL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
ERT1800	ERGOJECT™ 1.8-TL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses
ERT1600	ERGOJECT™ 1.6-TL Injector Set	1.6	1.8	2.0	1-piece lenses (hydrophilic only)

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.



ERGOJECT™

SINGLE-USE INJECTOR SETS, BACK-LOADED


Compared to classical syringe injectors, the pen-like operation of ERGOJECT™ provides a much higher stability since the operation point is around 2cm closer to the eye. Consequently, lever forces and coherent up- and downturns of the injector are significantly reduced.

Different from classical two-handed screw injectors, the second hand remains free and it can be used for stabilizing the eye or for manipulating the lens within the capsular bag.

Cartridges and top-loaded loading chambers are identical to the ones used in the ACCUJECT™ system. This makes the use of ERGOJECT™ very convenient for current ACCUJECT™ users.

BACK-LOADED / FOR HYDROPHOBIC LENSES

consisting of: ERGOJECT™ Injection system
with integrated cartridge and back-loaded loading chamber
(Box of 10)

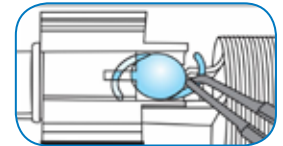

ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
ERB3000	ERGOJECT™ 3.0-BL Injector Set	2.4	2.7-2.8	3.1	1-piece lenses
ERB2600	ERGOJECT™ 2.6-BL Injector Set	2.2	2.5	2.8	1-piece lenses
ERB2400	ERGOJECT™ 2.4-BL Injector Set	2.1	2.4	2.7	1-piece lenses
ERB2200	ERGOJECT™ 2.2-BL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
ERB2100	ERGOJECT™ 2.1-BL Injector Set	1.9	2.1-2.2	2.4	1-piece lenses
ERB2000	ERGOJECT™ 2.0-BL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
ERB1800	ERGOJECT™ 1.8-BL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.



ACCUJECT™ DUAL



SINGLE-USE INJECTOR SETS, STERILE **NEW!**

With our single-use ACCUJECT™ DUAL injector, the surgeon is able to switch from either 100% push to 100% screw at any time during the IOL-injection procedure. Its movable finger flanges entirely eliminate usability issues. It is available as top-loading or back-loading version.

TOP-LOADED / FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: ACCUJECT™ DUAL Injection system
with integrated cartridge and top-loaded loading chamber
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
ADT3000	ACCUJECT™ DUAL 3.0-TL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
ADT2600	ACCUJECT™ DUAL 2.6-TL Injector Set	2.2	2.5	2.8	1-piece lenses
ADT2400	ACCUJECT™ DUAL 2.4-TL Injector Set	2.1	2.4	2.7	1-piece lenses
ADT2200	ACCUJECT™ DUAL 2.2-TL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
ADT2100	ACCUJECT™ DUAL 2.1-TL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
ADT2000	ACCUJECT™ DUAL 2.0-TL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
ADT1800	ACCUJECT™ DUAL 1.8-TL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses
ADT1600	ACCUJECT™ DUAL 1.6-TL Injector Set	1.6	1.8	2.0	1-piece lenses (hydrophilic only)

BACK-LOADED / FOR HYDROPHOBIC LENSES

consisting of: ACCUJECT™ DUAL Injection system
with integrated cartridge and back-loaded loading chamber
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
ADB3000	ACCUJECT™ DUAL 3.0-BL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
ADB2600	ACCUJECT™ DUAL 2.6-BL Injector Set	2.2	2.5	2.8	1-piece lenses
ADB2400	ACCUJECT™ DUAL 2.4-BL Injector Set	2.1	2.4	2.7	1-piece lenses
ADB2200	ACCUJECT™ DUAL 2.2-BL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
ADB2100	ACCUJECT™ DUAL 2.1-BL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
ADB2000	ACCUJECT™ DUAL 2.0-BL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
ADB1800	ACCUJECT™ DUAL 1.8-BL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.

ACCUJECT™

SINGLE-USE INJECTOR SETS, TOP-LOADED, STERILE

The compact ACCUJECT™ design with integrated cartridge enables a simple, predictable loading and positioning of any lens. It is designed for 1- and 3-piece lenses.

FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: ACCUJECT™ Injection system
with integrated cartridge and top-loaded loading chamber
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
LP604580	ACCUJECT™ 3.0-3P Injector Set	2.4	2.8	3.2	3-piece lenses
LP604520	ACCUJECT™ 2.5-3P Injector Set	2.2	2.5	2.8 - 2.9	3-piece lenses
LP604570	ACCUJECT™ 3.0-1P Injector Set	2.4	2.7-2.8	3.1	1-piece lenses
LP604500	ACCUJECT™ 2.6-1P Injector Set	2.2	2.5	2.8	1-piece lenses
LP604600	ACCUJECT™ 2.4-1P Injector Set	2.1	2.4	2.7	1-piece lenses
LP604530	ACCUJECT™ 2.2-1P Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
LP604590	ACCUJECT™ 2.1-1P Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
LP604510	ACCUJECT™ 2.0-1P Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
LP604540	ACCUJECT™ 1.8-1P Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses
LP604550	ACCUJECT™ 1.6-1P Injector Set	1.6	1.8	2.0	1-piece lenses (hydrophilic only)

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.



for 1-piece IOLs



for 3-piece IOLs

ACCUJECT™

SINGLE-USE INJECTOR SETS, BACK-LOADED

The intelligent ACCUJECT™ back-loaded technique combines the advantages of the silicone cushion technology with the advantages of the classical back-loaded injectors. It is therefore impossible to ride over or scratch the lens as often seen when using a hard-tip injector.

FOR HYDROPHOBIC LENSES

consisting of: ACCUJECT™ Injection system
with integrated cartridge and back-loaded loading chamber
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
LP604575	ACCUJECT™ 3.0-BL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
LP604505	ACCUJECT™ 2.6-BL Injector Set	2.2	2.5	2.8	1-piece lenses
LP604605	ACCUJECT™ 2.4-BL Injector Set	2.1	2.4	2.7	1-piece lenses
LP604535	ACCUJECT™ 2.2-BL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
LP604595	ACCUJECT™ 2.1-BL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
LP604515	ACCUJECT™ 2.0-BL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
LP604545	ACCUJECT™ 1.8-BL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.



ACCUJECT™ SCREW



SINGLE-USE INJECTOR SETS, STERILE

Our single-use ACCUJECT™ SCREW injection system allows a quick release of the lens from the loading chamber using the known push-mechanism. As soon as the screw thread is grasped, the lens can be injected smoothly and controlled via the screw advancement. It is available as top-loading or back-loading version.

TOP-LOADED / FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: ACCUJECT™ SCREW Injection system
with integrated cartridge and top-loaded loading chamber
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
AST3000	ACCUJECT™ SCREW 3.0-TL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
AST2400	ACCUJECT™ SCREW 2.4-TL Injector Set	2.1	2.4	2.7	1-piece lenses
AST2200	ACCUJECT™ SCREW 2.2-TL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
AST2100	ACCUJECT™ SCREW 2.1-TL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
AST2000	ACCUJECT™ SCREW 2.0-TL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
AST1800	ACCUJECT™ SCREW 1.8-TL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses (hydrophilic only)
AST1600	ACCUJECT™ SCREW 1.6-TL Injector Set	1.6	1.8	2.0	1-piece lenses (hydrophilic only)

BACK-LOADED / FOR HYDROPHOBIC LENSES

consisting of: ACCUJECT™ SCREW Injection system
with integrated cartridge and back-loaded loading chamber
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
ASB3000	ACCUJECT™ SCREW 3.0-BL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
ASB2600	ACCUJECT™ SCREW 2.6-BL Injector Set	2.2	2.5	2.8	1-piece lenses
ASB2400	ACCUJECT™ SCREW 2.4-BL Injector Set	2.1	2.4	2.7	1-piece lenses
ASB2200	ACCUJECT™ SCREW 2.2-BL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
ASB2100	ACCUJECT™ SCREW 2.1-BL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
ASB2000	ACCUJECT™ SCREW 2.0-BL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
ASB1800	ACCUJECT™ SCREW 1.8-BL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.

VISCOJECT™-BIO

SINGLE-USE INJECTOR SETS, STERILE

ENGINEERED BIOFILM coating makes VISCOJECT™-BIO the leading MICS lens delivery system a truly unique product. ENGINEERED BIOFILM offers outstanding gliding properties for IOLs to be safely injected through even smaller incision sizes.


FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: VISCOJECT™ Injector with silicone cushion
VISCOJECT™-BIO Cartridge
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
LP604335C	VISCOJECT™-BIO 2.7 Injector Set	2.2 - 2.3	2.5 - 2.6	2.9	1-piece lenses
LP604340C	VISCOJECT™-BIO 2.2 Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
LP604350C	VISCOJECT™-BIO 1.8 Injector Set	1.8	2.0	2.2	1-piece lenses (hydrophilic only)
LP604370C	VISCOJECT™-BIO 1.6 Injector Set	1.6	1.8	2.0	1-piece lenses (hydrophilic only)
LP604360C	VISCOJECT™-BIO 1.5 Injector Set	1.4	1.6	1.8	1-piece lenses (for micro-incision lenses only)

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.

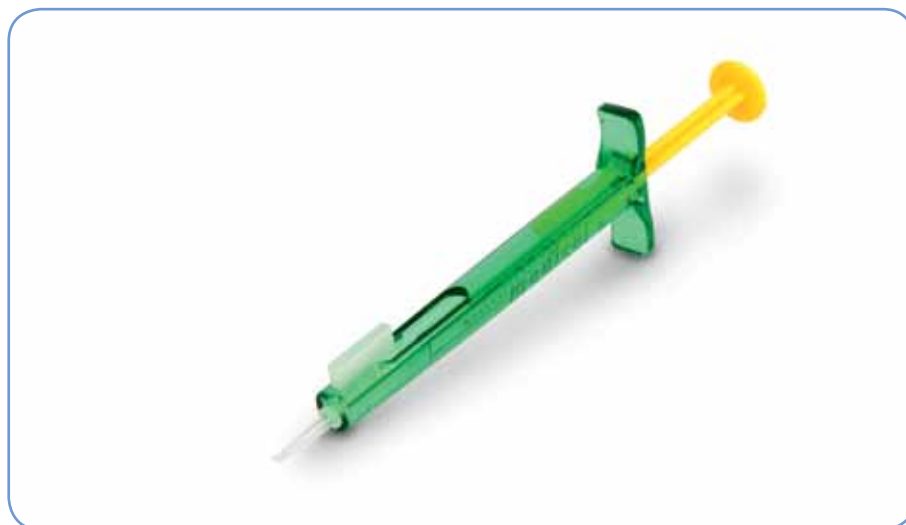

VISCOJECT™-BIO

SINGLE-USE CARTRIDGE SETS, STERILE

A wide range of cartridge sizes and its affordable pricing makes the VISCOJECT™-BIO with ENGINEERED BIOFILM coating the perfect allround system for single-piece IOLs of any design. Available in combination with a reusable titanium injector.

FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: VISCOJECT™-BIO Cartridge with silicone cushion
(Box of 10)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
LP604235C	VISCOJECT™-BIO 2.7 Cartridge Set	2.2 - 2.3	2.5 - 2.6	2.9	1-piece lenses
LP604240C	VISCOJECT™-BIO 2.2 Cartridge Set	1.9 - 2.0	2.2	2.5	1-piece lenses
LP604250C	VISCOJECT™-BIO 1.8 Cartridge Set	1.8	2.0	2.2	1-piece lenses (hydrophilic only)
LP604270C	VISCOJECT™-BIO 1.6 Cartridge Set	1.6	1.8	2.0	1-piece lenses (hydrophilic only)
LP604260C	VISCOJECT™-BIO 1.5 Cartridge Set	1.4	1.6	1.8	1-piece lenses (hydrophilic only)

*Given incision sizes are indications only based on average experiences.
The individually required incision sizes may vary depending on surgical techniques.

REUSABLE INJECTORS, NON-STERILE

ITEM CODE	PRODUCT NAME	TYPE	TO BE USED WITH VISCOJECT™-BIO CARTRIDGE SETS
LI604205	VISCOJECT™ Injector „ATRAUMATIC“	syringe type	LP604235C / LP604240C / LP604250C / LP604270C / LP604260C
LI604215	VISCOJECT™ Injector „ATRAUMATIC“	screw type	LP604235C / LP604240C / LP604250C / LP604270C / LP604260C

syringe type



screw type



VISCOJECT™ *easy*

SINGLE-USE INJECTOR SETS, STERILE

VISCOJECT™ *easy* is the youngest member of the VISCOJECT™ family. The loading of the lens is done with the aid of a loading unit, which defines the angular position of the cartridge wing and the management of the loading is considerably simplified. The design of the cartridge tip is equivalent to the advanced VISCOJECT™-BIO cartridge.

FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: VISCOJECT™ *easy* Injector with silicone cushion
 VISCOJECT™ *easy* Cartridge
 VISCOJECT™ *easy* Loading unit
 (Box of 10)



loading unit

ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
VY2900	VISCOJECT™ <i>easy</i> 2.9-1P Injector Set	2.4	2.7	3.0 - 3.1	1-piece lenses
VY2700	VISCOJECT™ <i>easy</i> 2.7-1P Injector Set	2.2 - 2.3	2.5 - 2.6	2.9	1-piece lenses
VY2400	VISCOJECT™ <i>easy</i> 2.4-1P Injector Set	2.1	2.4	2.7	1-piece lenses
VY2200	VISCOJECT™ <i>easy</i> 2.2-1P Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses

*Given incision sizes are indications only based on average experiences.
 Individual required incision sizes may vary depended on surgical techniques.



VISCOJECT™ *easy*

SINGLE-USE CARTRIDGE SETS, STERILE

VISCOJECT™ *easy* is also available in combination with reusable titanium injectors. The loading of the lens is done with the aid of a loading unit, which defines the angular position of the cartridge wing and the management of the loading is considerably simplified. The design of the cartridge tip is equivalent to the advanced VISCOJECT™-BIO cartridge.

FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: VISCOJECT™ *easy* Cartridge with silicone cushion
VISCOJECT™ *easy* Loading unit
(Box of 10)



loading unit

ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			LOADING UNIT YES / NO	IOL TYPE
		wound-assisted	into-the-wound	into-the-bag		
VC2900	VISCOJECT™ <i>easy</i> 2.9-1P Cartridge Set	2.4	2.7	3.0 - 3.1	Yes	1-piece lenses
VC2700	VISCOJECT™ <i>easy</i> 2.7-1P Cartridge Set	2.2 - 2.3	2.5 - 2.6	2.9	Yes	1-piece lenses
VC2400	VISCOJECT™ <i>easy</i> 2.4-1P Cartridge Set	2.1	2.4	2.7	Yes	1-piece lenses
VC2200	VISCOJECT™ <i>easy</i> 2.2-1P Cartridge Set	1.9 - 2.0	2.2	2.5	Yes	1-piece lenses

*Given incision sizes are indications only based on average experiences. Individual required incision sizes may vary depended on surgical techniques.

REUSABLE INJECTORS, NON-STERILE

ITEM CODE	PRODUCT NAME	TYPE	TO BE USED WITH VISCOJECT™ <i>easy</i> CARTRIDGE SETS
LI604305	VISCOJECT™ <i>easy</i> Injector „ATRAUMATIC“	syringe type	VC2900 / VC2700 / VC2400 / VC2200
LI604315	VISCOJECT™ <i>easy</i> Injector „ATRAUMATIC“	screw type	VC2900 / VC2700 / VC2400 / VC2200

syringe type



screw type



VISCOJECT™ **eco**

SINGLE-USE INJECTOR SETS, STERILE

Are you located in a price sensitive market? The attractive price performance ratio of the VISCOJECT™ **eco** is the perfect solution. Most importantly, equipped with the traditional VISCOJECT™ cartridge, the VISCOJECT™ **eco** does not compromise on quality and safety.

FOR HYDROPHOBIC AND HYDROPHILIC LENSES

consisting of: VISCOJECT™ **eco** Injector with silicone cushion
VISCOJECT™ Cartridge
(Box of 20)



ITEM CODE	PRODUCT NAME	INCISION SIZE* (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
VE2700	VISCOJECT™ eco 2.7 Injector Set	2.3	2.6 - 2.7	2.9 - 3.0	1-piece lenses
VE2200	VISCOJECT™ eco 2.2 Injector Set	2.1 - 2.2	2.4	2.7 - 2.8	1-piece lenses

*Given incision sizes are indications only based on average experiences.
Individual required incision sizes may vary depended on surgical techniques



MULTIJECT™

SINGLE-USE CARTRIDGES, STERILE

The combination of the MULTIJECT™ reusable titanium lens injector and the MULTIJECT™ cartridges results in a sound and well proven complete lens injection system compatible with silicone as well as acrylic lenses of any design. The MULTIJECT™ allows injection of any style of foldable 1- and 3-piece acrylic and silicone lenses.

FOR HYDROPHOBIC AND HYDROPHILIC LENSES


consisting of: MULTIJECT™ Cartridge



ITEM CODE	PRODUCT NAME	IOL TYPE	INCISION SIZE*
LC604215	MACROGLIDE™ Cartridge	1- and 3-piece lenses	sub 3.5mm
LC604210	OVALGLIDE™ Cartridge	1- and 3-piece lenses	sub 3.2mm
LC604220M	MICROGLIDE™ Cartridge	1- and 3-piece lenses	sub 3.0mm
LS604500	SILICONE CUSHION	to be used with reusable injector LI604205 / LI604215	

*Given incision sizes are indications only based on average experiences. Individual required incision sizes may vary depended on surgical techniques

REUSABLE INJECTORS, NON-STERILE

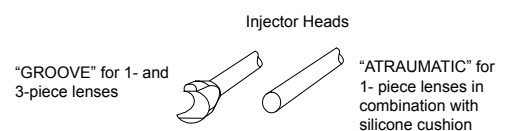

ITEM CODE	PRODUCT NAME	TYPE	IOL TYPE	TO BE USED WITH MULTIJECT™ CARTRIDGES
LI604200	MULTIJECT™ Injector „GROOVE“	syringe type	1- and 3-piece lenses	LC604215 LC604210 LC604220M
LI604205	MULTIJECT™ Injector „ATRAUMATIC“	syringe type	1-piece lenses	
LI604210	MULTIJECT™ Injector „GROOVE“	screw type	1- and 3-piece lenses	
LI604215	MULTIJECT™ Injector „ATRAUMATIC“	screw type	1-piece lenses	

syringe type



screw type



I/A INSTRUMENTS

SINGLE-USE SWISS COAXIAL I/A SYSTEMS, STERILE

Ergonomic grip design!

Beside the Classic Design – known and valued by surgeons all-over the world – Medicel offers a complete new range of disposable I/A handles. The new I/A handles. The new design responds to surgeons demand to improve ergonomics and grip of I/A handles without compromising product quality. Latest manufacturing technology allows us to offer an attractive priced product without compromising the quality.

consisting of: Coaxial I/A Instrument
(Box of 10)



ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD
SMS325P	20G / 0.9mm	2.5mm	angled 45°	0.35mm	polished
SMS335	20G / 0.9mm	2.5mm	angled 65°	0.35mm	textured
SMS355P	21G / 0.8mm	2.2mm	angled 45°	0.35mm	polished
SMS370P	23G / 0.65mm	1.8mm	angled 45°	0.30mm	polished
SMS410 (FLACS)	20G / 0.9mm	2.5mm	30° bevelled tip	0.60mm	n/a
SMS430 (FLACS)	21G / 0.8mm	2.2mm	30° bevelled tip	0.60mm	n/a

WITH SILICONE TIP



ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD
SMS390S	21G / 0.8mm	2.2mm	angled 45°	0.40mm	silicone tip

SINGLE-USE SWISS BIMANUAL I/A SYSTEMS, STERILE

consisting of: Bimanual I/A Instrument
(Box of 10)



ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD	IRRIGATION CANNULA	IRRIGATION PORT
SBS310	21G / 0.8mm	1.1mm	bullet shape	0.35mm	textured	bullet shape	2 ports, 0.45mm
SBS320	23G / 0.65mm	0.8mm	bullet shape	0.30mm	textured	bullet shape	2 ports, 0.35mm

I/A INSTRUMENTS

SINGLE-USE SWISS COAXIAL I/A SYSTEMS, STERILE


The Medicel range of coaxial irrigation / aspiration systems features straight and angled versions at 45° and 65°. All models are available as 19G, 20G, 21G and 23G version required for the coaxial micro-incision method. The textured aspiration tip can be used for selective polishing of the posterior capsule.

As all I/A systems are equipped with perfectly matching infusion sleeves, dimensional conditions are identical to the previous phaco-procedure and thus ensuring a smooth and totally leak-free cortical removal.

consisting of: Coaxial I/A Instrument
(Box of 10)



ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD
SMS120	19G / 1.1mm	2.8mm	angled 45°	0.35mm	textured
SMS130	19G / 1.1mm	2.8mm	angled 65°	0.35mm	textured
SMS115	20G / 0.9mm	2.5mm	straight	0.35mm	textured
SMS125P	20G / 0.9mm	2.5mm	angled 45°	0.35mm	polished
SMS135	20G / 0.9mm	2.5mm	angled 65°	0.35mm	textured
SMS155P	21G / 0.8mm	2.2mm	angled 45°	0.35mm	polished
SMS160P	21G / 0.8mm	2.2mm	curved	0.35mm	polished
SMS170P	23G / 0.65mm	1.8mm	angled 45°	0.30mm	polished
SMS175P	23G / 0.65mm	1.8mm	curved	0.30mm	polished
SMS180P	25G / 0.6mm	1.5mm	curved	0.30mm	polished

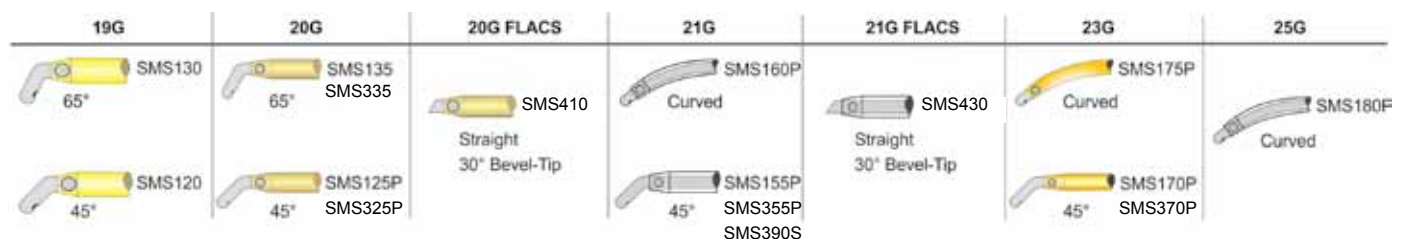

I/A INSTRUMENTS

SINGLE-USE SWISS BIMANUAL I/A SYSTEMS, STERILE

The curved irrigation and aspiration cannulas with elegant aeroplane-nose design have a diameter of just 0.8 mm (21G) and can be introduced effortlessly through the smallest paracentesis. At the same time, the inner lumens are designed with maximum size for efficient capsule cleaning with a stable anterior chamber.

The tip of the aspiration cannula is textured for capsule polishing. On the SBS105 the irrigation cannula is slightly pointed and open at the front. This ensures an active rinsing function.

consisting of: Bimanual I/A Instrument
(Box of 10)



ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD	IRRIGATION CANNULA	IRRIGATION PORT
SBS105	21G / 0.8mm	1.1mm	curved, aeroplane shape	0.35mm	textured	open	1 port, 0.5mm
SBS110A	21G / 0.8mm	1.1mm	curved, aeroplane shape	0.35mm	textured	curved, aeroplane shape	2 ports, 0.45mm
SBS122A	22G / 0.7mm	1.0mm	aeroplane shape	0.35mm	textured	curved, aeroplane shape	2 ports, 0.45mm
SBS123A	23G / 0.65mm	0.8mm	aeroplane shape	0.30mm	textured	curved, aeroplane shape	2 ports, 0.35mm

MULTI-USE BIMANUAL I/A SYSTEMS, STERILE (AUTOCLAVABLE UP TO 10 TIMES)

This unique semi-disposable line is fully autoclavable and validated for up to 10 sterilisations. Same as the single-use line, the multiple-use products are sterile packed.


Sets delivered sterile, components can be autoclaved at 121°C/20 min or 134°C/5 min up to 10 times.

consisting of: Bimanual I/A Instrument
(Box of 10)




ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD	IRRIGATION CANNULA	IRRIGATION PORT
SBS105RU	21G / 0.8mm	1.1mm	curved, aeroplane shape	0.35mm	textured	open	1 port, 0.5mm
SBS110ARU	21G / 0.8mm	1.1mm	curved, aeroplane shape	0.35mm	textured	curved, aeroplane shape	2 ports, 0.45mm

SBS105
SBS105RU

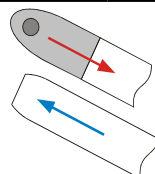

SBS110A
SBS110ARU

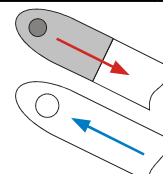

I/A INSTRUMENTS

REUSABLE SWISS BIMANUAL I/A SYSTEMS, NON-STERILE

The Medicel bimanual I/A systems can be combined in various configurations. All models are made of titanium. Thanks to the aeroplane-nose design, cannulas can be introduced easily through the smallest of incisions.

consisting of: Aspiration Instrument
(Box of 1)

I/A ASPIRATION INSTRUMENTS



ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD
SBA100A	21G / 0.8mm	1.1mm	aeroplane shape	0.35mm	polished
SBA103A	21G / 0.8mm	1.1mm	aeroplane shape	0.40mm	polished
SBA104A	21G / 0.8mm	1.1mm	aeroplane shape	0.40mm	textured
SBA105A	21G / 0.8mm	1.1mm	aeroplane shape	0.35mm	textured
SBA101B	23G / 0.6mm	0.8mm	bullet shape	0.30mm	polished
SBA106B	23G / 0.6mm	0.8mm	bullet shape	0.30mm	textured

consisting of: Irrigation Instrument
(Box of 1)

I/A IRRIGATION INSTRUMENTS



ITEM CODE	SIZE	INCISION SIZE	IRRIGATION CANNULA	IRRIGATION PORT	IRRIGATION HEAD
SBI120A	21G / 0.8mm	1.1mm	aeroplane shape	2 ports, 0.45mm	polished
SBI125A	21G / 0.8mm	1.1mm	aeroplane shape	2 ports, 0.45mm	textured
SBI121B	23G / 0.6mm	0.8mm	bullet shape	2 ports, 0.30mm	polished
SBI126B	23G / 0.6mm	0.8mm	bullet shape	2 ports, 0.30mm	textured

Bullet-shape



aeroplane-shape



4K F-SONIC PHACO SETS

4K F-SONIC FLARED PHACO ACCESSORIES SETS



Stable chamber even in high vacuum levels, efficient emulsification to reduce effective phaco-time without compromising safety - this is our latest member of the phaco tip family - F-SONIC.

compatible with the following phaco machines: Bausch&Lomb, AMO/J&J, Nidek, Oertli, Ruck/Croma, Geuder, DORC

Single-use sets (sterile) consisting of 10 titanium tips with 10 tip-wrenches / 10 infusion sleeves / 10 test chambers

ITEM CODE	PRODUCT NAME	SIZE	AVAILABLE DEGREE	INCISION SIZE
PA603019FK	4K F-SONIC 19G/30°	19G / 1.1mm	30° (Standard)	2.8mm
PA603020FK	4K F-SONIC 20G/30°	20G / 0.9mm	30° (Standard)	2.5mm
PA603021FK	4K F-SONIC 21G/30°	21G / 0.8mm	30° (Standard)	2.2mm
PA603023FK	4K F-SONIC 23G/30°	23G / 0.7mm	30° (Standard)	1.8mm

10 times multi-use (sterile) sets consisting of 10 titanium tips with 10 autoclavable tip-wrenches and infusion sleeves

ITEM CODE	PRODUCT NAME	SIZE	AVAILABLE DEGREE	INCISION SIZE
PA603019FKRU	4K F-SONIC 19G/30°	19G / 1.1mm	30° (Standard)	2.8mm
PA603020FKRU	4K F-SONIC 20G/30°	20G / 0.9mm	30° (Standard)	2.5mm
PA603021FKRU	4K F-SONIC 21G/30°	21G / 0.8mm	30° (Standard)	2.2mm
PA603023FKRU	4K F-SONIC 23G/30°	23G / 0.7mm	30° (Standard)	1.8mm



4K PHACO ACCESSORIES

4K PHACO ACCESSORIES SETS



Like the I/A systems, phaco accessories are available as single-use or reusable versions. The full range of phaco tips is compatible with various systems and is made of high-grade titanium. Medical-grade silicone sleeves ensure best cooling of the phaco tip, protection of incision and stable anterior chamber during surgery, even when using high vacuum levels.

compatible with the following phaco machines: Bausch&Lomb, AMO/J&J, Nidek, Oertli, Ruck/Croma, Geuder, DORC

PRODUCT NAME	SIZE	AVAILABLE DEGREE	INCISION SIZE
CLASSIC	19G / 1.1mm	30° (Standard)*	2.8mm
MICROSONIC	20G / 0.9mm	30° (Standard)*	2.5mm
HYPERSONIC	19/20G / 1.1mm	30° (Standard)*	2.8mm
NANO	20/21G / 0.8mm	30° (Standard)*	2.2mm
SUBSONIC	21G / 0.8mm	30° (Standard)*	2.2mm
COMICS-SUB2	23G / 0.7mm	30° (Standard)*	1.8mm

*degrees 45° / 15° / 0° on demand

Single-use sets (sterile) consisting of 10 titanium tips with 10 tip-wrenches / 10 infusion sleeves / 10 test chambers

CLASSIC 19G, 30°	PA603500	
MICROSONIC 20G, 30°	PA603510	
HYPERSONIC 19/20G, 30°	PA603520	
NANO 20/21G, 30° (flared)	PA603530	
SUBSONIC 21G, 30°	PA603540	
COMICS-SUB2 23G, 30°	PA603550	

10 times multi-use (sterile) sets consisting of 10 titanium tips with 10 autoclavable tip-wrenches and infusion sleeves

CLASSIC 19G, 30°	PA603500RU	
MICROSONIC 20G, 30°	PA603510RU	
HYPERSONIC 19/20G, 30°	PA603520RU	
NANO 20/21G, 30° (flared)	PA603530RU	
SUBSONIC 21G, 30°	PA603540RU	
COMICS-SUB2 23G, 30°	PA603550RU	

Reusable phaco tips (non-sterile) consisting of titanium 6 tips

CLASSIC 19G, 30°	PT602500B	
MICROSONIC 20G, 30°	PT602510B	
HYPERSONIC 19/20G, 30°	PT602520B	
NANO 20/21G, 30° (flared)	PT602530B	
SUBSONIC 21G, 30°	PT602540B	
COMICS-SUB2 23G, 30°	PT602550B	
PHACO TIP WRENCH (4K)	PK602401	

5R PHACO ACCESSORIES

5R PHACO ACCESSORIES SETS



Like the I/A systems, phaco accessories are available as single-use or reusable versions. The full range of phaco tips is compatible with various systems and is made of high-grade titanium. Medical-grade silicone sleeves ensure best cooling of the phaco tip, protection of incision and stable anterior chamber during surgery, even when using high vacuum levels.

compatible with the following phaco machines: Bausch&Lomb, AMO/J&J, Nidek, Oertli, Ruck/Croma, Geuder, DORC

PRODUCT NAME	SIZE	AVAILABLE DEGREE	INCISION SIZE
CLASSIC	19G / 1.1mm	30° (Standard)*	2.8mm
MICROSONIC	20G / 0.9mm	30° (Standard)*	2.5mm
HYPERSONIC	19/20G / 1.1mm	30° (Standard)*	2.8mm
NANO	20/21G / 0.8mm	30° (Standard)*	2.2mm
SUBSONIC	21G / 0.8mm	30° (Standard)*	2.2mm
COMICS-SUB2	23G / 0.7mm	30° (Standard)*	1.8mm

*degrees 45° / 15° / 0° on demand

Single-use sets (sterile) consisting of 10 titanium tips with 10 tip-wrenches / 10 infusion sleeves / 10 test chambers

CLASSIC 19G, 30°	PA603000	
MICROSONIC 20G, 30°	PA603100	
HYPERSONIC 19/20G, 30°	PA603200	
NANO 20/21G, 30° (flared)	PA603190	
SUBSONIC 21G, 30°	PA603070	
COMICS-SUB2 23G, 30°	PA603020	

10 times multi-use (sterile) sets consisting of 10 titanium tips with 10 autoclavable tip-wrenches and infusion sleeves

CLASSIC 19G, 30°	PA603000RU	
MICROSONIC 20G, 30°	PA603100RU	
HYPERSONIC 19/20G, 30°	PA603200RU	
NANO 20/21G, 30° (flared)	PA603190RU	
SUBSONIC 21G, 30°	PA603070RU	
COMICS-SUB2 23G, 30°	PA603020RU	

Reusable phaco tips (non-sterile) consisting of titanium 6 tips

CLASSIC 19G, 30°	PT602200B	
MICROSONIC 20G, 30°	PT602210B	
HYPERSONIC 19/20G, 30°	PT602250B	
NANO 20/21G, 30° (flared)	PT602240B	
SUBSONIC 21G, 30°	PT602070B	
COMICS-SUB2 23G, 30°	PT602220B	
PHACO TIP WRENCH (5R)	PK602402	

XXK PHACO ACCESSORIES

XXK PHACO ACCESSORIES SETS



SUBSONIC Phaco tip sets are now available in the new version suitable for Alcon Infinity® in 21G straight and curved version.

compatible with the following phaco machines: Alcon Infinity®, Constellation®, Centurion®

Single-use sets (sterile) consisting of 10 titanium tips with 10 tip-wrenches / 10 infusion sleeves / 10 test chambers


ITEM CODE	PRODUCT NAME	SIZE	AVAILABLE DEGREE	INCISION SIZE
PA603391	SUBSONIC 21G	21G / 0.8mm	30° (Standard), straight	2.2mm
PA603391C	SUBSONIC 21G	21G / 0.8mm	30° (Standard), curved	2.2mm



PHACO SLEEVES

INFUSION SLEEVES AND TEST CHAMBER

Medicel COOLFLOW infusion sleeves are made from medical-grade silicone and are manufactured on high-precision moulds in-house. COOLFLOW infusion sleeves are designed to ensure maximal thermal protection to the cornea whilst maintaining a close and stable chamber.

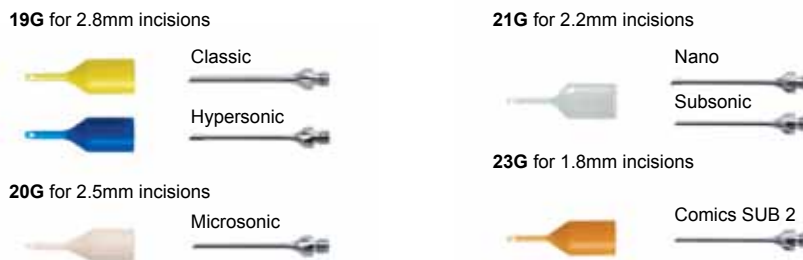

Single-use sets (sterile) consisting of 10 infusion sleeves / 10 test chambers

ITEM CODE	PRODUCT NAME	COLOUR	SIZE	INCISION SIZE	TO BE USED WITH THE FOLLOWING PHACO TIPS
PA603901	ACCESSORIES SET 19G	BLUE	19G / 1.1mm	2.8mm	CLASSIC / HYPERSONIC
PA603902	ACCESSORIES SET 20G	WHITE	20G / 0.9mm	2.5mm	MICROSONIC
PA603903	ACCESSORIES SET 21G	CLEAR	21G / 0.8mm	2.2mm	NANO / SUBSONIC
PA603904	ACCESSORIES SET 23G	ORANGE	23G / 0.7mm	1.8mm	COMICS SUB-2

Reusable phaco infusion sleeves (non-sterile) consisting of 6 infusion sleeves

ITEM CODE	PRODUCT NAME	COLOUR	SIZE	INCISION SIZE	TO BE USED WITH THE FOLLOWING PHACO TIPS
CF602300	COOLFLOW 19G	YELLOW	19G / 1.1mm	2.8mm	CLASSIC / HYPERSONIC
CF602301	COOLFLOW 19G	BLUE	19G / 1.1mm	2.8mm	CLASSIC / HYPERSONIC
CF602310	COOLFLOW 20G	WHITE	20G / 0.9mm	2.5mm	MICROSONIC
CF602340	COOLFLOW 21G	CLEAR	21G / 0.8mm	2.2mm	SUBSONIC
CF602350	COOLFLOW 23G	ORANGE	23G / 0.8mm	2.0mm	NANO
TC602450	TEST CHAMBER	TRANSPARENT	n/a	n/a	n/a

Reusable phaco infusion sleeves (non-sterile) consisting of 4 infusion sleeves

ITEM CODE	PRODUCT NAME	COLOUR	SIZE	INCISION SIZE	TO BE USED WITH THE FOLLOWING PHACO TIPS
CF602320	COOLFLOW 19G, rigid no-burn	YELLOW	19G / 1.1mm	2.8mm	CLASSIC / HYPERSONIC
CF602330	COOLFLOW 20G rigid no-burn	WHITE	20G / 1.1mm	2.5mm	MICROSONIC

COMPATIBLES

SINGLE-USE ULTRA-HIGHVAC XXK CASSETTE SYSTEMS, STERILE

consisting of: 1 XXK ULTRA-HIGHVAC cassette
1 I/A tubing set
1 Administration set (with filter)
1 Tray support cover
(Box of 12)



compatible with the following phaco machines: Alcon Legacy® Series 20000

ITEM CODE	CASSETTE SYSTEM	PHACO TIPS	PHACO SLEEVES	TEST CHAMBER	INCISION SIZE
XXK100	XXK ULTRA-HIGHVAC set	not included	not included	not included	n/a
XXK101	XXK ULTRA-HIGHVAC set	not included	19G sleeve, blue	included	2.8mm
XXK120	XXK ULTRA-HIGHVAC set	MICROSONIC 20G, 30°	20G sleeve, white	included	2.5mm

DMEK INJECTORS

SINGLE-USE DMEK INJECTOR SETS, STERILE

For safe implantation of grafts during Descemet Membrane Endothelial Keratoplasty

Medicel introduces the first CE approved injector specifically designed for the safe implantation of grafts during Descemet Membrane Endothelial Keratoplasty.

consisting of: ENDOJECT™ Injector with silicone cushion
ENDOJECT™ Cartridge
(Box of 1)



ITEM CODE	PRODUCT NAME	INCISION SIZE	SILICONE CUSHION
EJ2200	ENDOJECT™ 2.2	2.5mm	Yes

According to Prof. Dr. Friedrich Kruse

consisting of: GÜELL DMEK Injector
with integrated cartridge
(Box of 1)



ITEM CODE	PRODUCT NAME	INCISION SIZE	SILICONE CUSHION
DK2200	GÜELL DMEK 2.2	2.2mm	Yes

According to Dr. José Güell