


COAXIAL SOFT-TIP I/A SYSTEM

For ultra safe irrigation and aspiration
during cataract surgery

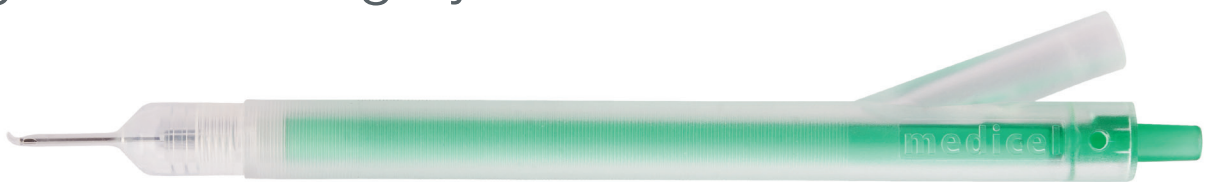

Distributed by:
Medix21 Surgical
29 Jarden Mile Ngauranga 6035
p 0800 633 4921 / 021 734 762
e camille@medix21.co.nz
w medix21surgical.co.nz

Mediceal AG Tel. +41 71 727 10 50
Dornierstrasse 11 Fax +41 71 727 10 55
CH-9423 Altenrhein info@mediceal.com
SWITZERLAND www.mediceal.com


mediceal
SWISS TECHNOLOGY FOR SURGERY

COAXIAL SOFT-TIP I/A SYSTEM

The Soft-Tip I/A enhances safety during cataract surgery by including two elements: Visibility and smooth touch.

- Touching the thin capsular bag with the metal tip of conventional instruments generates neither haptic nor optic feedback to the surgeon. The silicone tip of Soft-Tip I/A can bend once pushed to the capsular bag. The visible bending tip indicates the contact with the capsular bag to the surgeon.
- The adhesive material properties of silicone ease the capsular polishing procedure and reduce the risk of capsular rupture.

Single-use I/A system

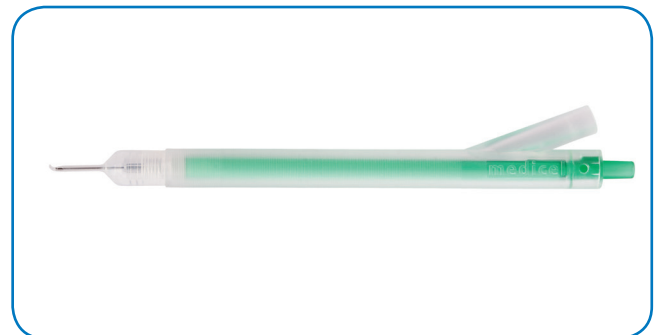

- Soft-tip
- Reduced risk of capsular rupture
- Safe lens manipulation

Single-use sets (sterile)

ITEM CODE	SIZE	INCISION SIZE	ASPIRATION CANNULA	ASPIRATION PORT	ASPIRATION HEAD
SMS390S	21G / 0.8mm	2.2mm	angled 45°	0.40mm	silicone tip